

BM3VET User Manual

Code 33719 Veterinary Monitor
Rev. 2.4



Rev. 2.4

Table of Contents

BM3VET User Manual	0
Table of Contents	1
1. BASIC	7
1.1 CE Standard Information	8
1.2 Read before Use	9
How to Contact Us	Errore. Il segnalibro non è definito.
Warranty Period	10
Warning, Caution, Note.....	11
General Precaution on Environment.....	12
General Precaution on Electric Safety	17
Cleaning Applied Parts	19
1.3 Product Components	21
Product Outline	21
Principal Characters of Product	21
Product Configuration	22
Option Product	22
Product Body Configuration	23
Accessories	25
Equipment Sign.....	26
1.4 Function and Key	29
External Function	29
Operation Key	29
1.5 Standard Power Supply Application	31
DC Power.....	31
1.6 Battery Power Supply Application	32
Operation	32
The Impact of Lithium-Ion Battery Technology on the Battery	34
Conditioning Guideline.....	34

Storage Guideline	34
How to Recycle the Battery.....	34
1.7 General Manu Operation	35
Screen Composition.....	35
Menu Selection	36
Menu Composition	36
2. ANIMAL/DATA MANAGEMENT.....	40
2.1 ADMIT.....	42
ANIMAL TYPE	42
CHANGE ANIMAL INFORMATION.....	43
DEFAULTS SETTING.....	43
HEIGHT.....	44
WEIGHT.....	44
2.2 ALARM.....	45
Alarm for the Product.....	46
ALL LIMITS	47
ALARM PRINT.....	47
ALARM VOLUME.....	48
ALARM LEVEL.....	48
PARAMETER LEVEL	49
ARRHYTHM LEVEL.....	49
ALARM REVIEW.....	50
ALARM LIST	51
SAVING CONDITION	52
NURSE CALL.....	53
3. SETUP.....	54
3.1 SETUP	55
DISPLAY	55
SET PARA	56
WAVE SELECT	56
SET DATE & TIME.....	57
SET TIME.....	57
SET DATE	58

HR SOURCE.....	58
SWEEP SPEED.....	59
KEY SOUND.....	59
DEMO.....	59
USER SERVICE.....	60
SET UNIT NAME.....	60
SET BED NUMBER.....	61
AC FILTER.....	61
SYSTEM.....	62
MAKER SERVICE.....	62
4. TREND.....	63
4.1 TREND.....	64
GRAPHIC TREND.....	65
TIME PERIOD.....	66
TABULAR TREND.....	67
TIME INTERVAL.....	68
TREND WINDOW SETUP.....	68
TIME PERIOD.....	69
SET TREND PARA.....	70
TREND PRINT.....	70
5. ECG.....	71
5.1 Introduction.....	72
Colors and Standards of Cables.....	72
Position of ECG Connector and Measuring Cable.....	72
Attaching Electrodes to the Animal.....	75
Choosing an ECG lead for Arrhythmia Monitoring.....	76
Information on the ECG waveform.....	76
Position of 3-Lead Wire Electrode.....	77
Position of 5-Lead Wire Electrode.....	78
5.2 ECG Data Window.....	79
5.3 ECG Data Setup.....	82
TRACE 1 LEAD SELECT.....	82
ALARM LIMIT.....	83

ALARM SOUND	84
QRS VOLUME	85
DISPLAY	85
ECG SPEED	86
ECG SIZE	86
HR SOURCE.....	87
ANALYSIS SETTING	87
6. SpO₂.....	98
6.1 Outline.....	99
SpO2 Connector Location and Measuring Cable	99
6.2 SpO2 Data Window	101
Signal and Data Validity.....	102
6.3 SpO₂ Data Setup.....	104
RATE VOLUME	104
ALARM.....	105
ALARM LIMIT.....	105
ALARM SOUND.....	106
LEAD FAULT Condition	106
SPO2 Messages	107
7. RESPIRATION.....	108
7.1 Outline.....	109
7.2 Respiration Data Window.....	112
7.3 Respiration Data Setup.....	113
RESPIRATION SPEED	113
RESPIRATION.....	114
APNEA DETECT.....	114
ALARM.....	115
ALARM LIMIT.....	115
ALARM SOUND.....	116
8. NIBP.....	117
8.1 Outline.....	118
8.2 NIBP Data Window	121

8.3 NIBP Data Setup	122
ALARM.....	122
ALARM LIMIT.....	123
ALARM SOUND.....	124
CUFF SIZE.....	124
UNIT SELECT.....	125
INTERVAL.....	125
INFLATION	126
9. TEMPERATURE	127
9.1 Outline	128
9.2 Temperature Data Window	129
9.3 Temperature Data Setup	130
ALARM.....	130
ALARM LIMIT.....	130
ALARM SOUND.....	131
UNIT SELECT.....	132
10. PRINT	133
10.1 Print	134
Printer and Heat Sensitivity Paper.....	134
Function and Setup Menu.....	135
10.2 Paper Change	138
11. MESSAGE LIST	139
12. DEFAULT SETTING VALUE	140
1. HORSE Mode	140
2. DOG Mode	142
3. PUPPY Mode	144
4. CAT Mode	146
13. TROUBLE SHOOTING	148
1. Noise in ECG	148
2. SpO2 malfunction	149

3. Temp malfunction	149
4. NIBP malfunction	150
5. Abnormality in NIBP measurements	150
6. Failure in battery recharge	151
7. Power failure.....	152
8. Periodic noises.....	153
9. Print failure	154
14. SPECIFICATION	155
Ease of use	156
Additional Function	156
Monitor Environmental Specifications	156
Power	156
Monitor Performance Specifications	156
Graphical and Tabular Trends	157
ECG capacity	157
SpO ₂ capacity.....	157
Respiration Performance Specifications	158
NIBP capacity.....	158
Temperature Unit Performance Specifications.....	158
Accessories Included:	158
Option.....	159
Abbreviations and Symbols.....	160
PRODUCT WARRANTY	<i>Errore. Il segnalibro non è definito.</i>

1. BASIC

1.1 CE Standard Information

1.2 Read before Use

Warranty Period

Warning, Caution, Note

General Precaution on Environment

General Precaution on Electric Safety

Equipment Connection, Maintenance & Washing Equipment Connection

1.3 Product Components

Product Outline

Principal Characteristics of Product

Product Configuration and Option Product

Product Body Configuration

1.4 Function and Key

External Function

Operation Key

1.5 Standard Power Supply Application

1.6 Battery Power Supply Application

1.7 General Menu Operation

Screen Composition

Menu Selection

Menu Composition

1.1 CE Standard Information

Electromechanical safety standards met:

- EN 60601-1: 1990 + A1:1993 + A2: 1995 Medical Electrical Equipment, Part 1, General Requirements for Safety.
- IEC/EN 60601-1-2 :2001 Electromagnetic compatibility -Requirements and tests.
- EN 1060-1:1995 Non-invasive sphygmomanometers - Part 1: General requirements
- EN 1060-3:1997 Non-invasive sphygmomanometers - Part 3: Supplementary requirements for electro-mechanical blood pressure measuring systems
- EN ISO 9919:2005 Medical electrical equipment - Particular requirements for the basic safety and essential performance of pulse oximeter equipment for medical use (ISO 9919:2005)
- EN 60601-2-27:2006 Medical electrical equipment - Part 2-27: Particular requirements for the safety, including essential performance, of electrocardiographic monitoring equipment
- EN 60601-2-30:2000 Medical electrical equipment - Part 2-30: Particular requirements for the safety, including essential performance, of automatic cycling non-invasive blood pressure monitoring equipment
- EN 12470-4:2000 Clinical thermometers - Part 4: Performance of electrical thermometers for continuous measurement
- EN 60601-2-49:2001 Medical electrical equipment - Part 2-49: Particular requirements for the safety of multifunction animal monitoring equipment

1.2 Read before Use

GIMA services are always available to you.

In the event of malfunction or failure, contact us along with the model name, serial number, and product name of the equipment.

Warranty Period

- This product is manufactured and passed through strict quality control and thorough inspection.
- We provide a 1-year warranty period.
- We will repair or replace any part of the BM3Vet found to be defective in usual operating circumstance at no cost to you.
- This warranty does not apply to any defect caused by improper use, abuse, misuse or improper handling.

Warning, Caution, Note

For special emphasis on agreement, terms are defined as listed below in user manual. Users should operate the equipment according to all the warnings and cautions.



Warning
To inform that it may cause serious injury or death to the animal, property damage, material losses

Caution
To inform that it may cause no harm in life but lead to injury

Note
To inform that it is not dangerous but important “note” sign for proper installation, operation, and maintenance of the equipment.

General Precaution on Environment

- Do not keep or operate the equipment in the environment listed below.

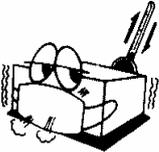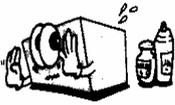	<p>Avoid placing in an area exposed to moisture. Do not touch the equipment with wet hands.</p>		<p>Avoid exposure to direct sunlight</p>
	<p>Avoid placing in an area where there is a high variation of temperature. Operating temperature ranges from 10(C to 40(C. Operating humidity ranges from 30% to 85%.</p>		<p>Do not use or store in the vicinity of an Electric heater</p>
	<p>Avoid placing in an area where there is an excessive humidity or poor ventilation.</p>		<p>Avoid placing in an area where there is an excessive shocks or vibration.</p>
	<p>Avoid placing in an area where chemicals are stored or where there is danger of gas leakage.</p>		<p>Avoid inserting dust and especially metal material into the equipment</p>
	<p>Do not open or disassemble the equipment. Bionet accepts no responsibility for unauthorized tampering service or repair.</p>		<p>Power off when the equipment is not fully installed. Otherwise, equipment could be damaged.</p>

CAUTIONS

Before Installation

Compatibility is critical to safe and effective use of this device. Please contact your local sales or service representative prior to installation to verify equipment compatibility.

Defibrillator Precaution

Animal signal inputs labeled with the CF and BF symbols with paddles are protected against damage resulting from defibrillation voltages. To ensure proper defibrillator protection, use only the recommended cables and lead wires.

Proper placement of defibrillator paddles in relation to the electrodes is required to ensure successful defibrillation.

Disposables

Disposable devices are intended for single use only. They should not be reused as performance could degrade or contamination could occur.

Disposal of your old appliance



1. When this crossed out wheeled bin symbol is attached to a product it means the product is covered by the European Directive 2002/96/EC.
2. All electrical and electronic products should be disposed of separately from the municipal waste stream via designated collection facilities appointed by the government or the local authorities.
3. The correct disposal of your old appliance will help prevent potential negative consequences for the environment and human health.
4. For more detailed information about disposal of your old appliance, please contact your city office, waste disposal service or the shop where you purchased the product.

Electrocute Precautions

To prevent skin burns, apply electrocute electrodes as far as possible from all other electrodes, a distance of at least 15 cm/6 in. is recommended.

EMC

Magnetic and electrical fields are capable of interfering with the proper performance of the device.

For this reason make sure that all external devices operated in the vicinity of the monitor comply with the relevant EMC requirements. X-ray equipment or MRI devices are possible sources of interference as they may emit higher levels of electromagnetic radiation.

Also, keep cellular phones and other telecommunication equipment away from the monitor.

CAUTIONS

Instruction for Use

For continued safe use of this equipment, it is necessary that the instructions are followed. However, instructions listed in this in no way supersede established medical practices concerning animal care.

Loss of Data

Should the monitor at any time temporarily lose animal data, the potential exists that active monitoring is not being done. Close animal observation or alternate monitoring devices should be used until monitor function is restored.

If the monitor does not automatically resume operation within 60 seconds, power cycle the monitor using the power on/off switch. Once monitoring is restored, you should verify correct monitoring state and alarm function.

Maintenance

Regular preventive maintenance should be carried out annually (Technical inspections). You are responsible for any requirements specific to your country.

MPSO

The use of a multiple portable socket outlet (MPSO) for a system will result in an enclosure leakage current equal to the sum of all individual earth leakage currents of the system if there is an interruption of the MPSO protective earth conductor. Do not use an additional extension cable with the MPSO as it will increase the chance of the single protective earth conductor interruption.

Negligence

BIONET does not assume responsibility for damage to the equipment caused by improperly vented cabinets, improper or faulty power, or insufficient wall strength to support equipment mounted on such walls.

NOTES

Power Requirements

Before connecting the device to the power line, check that the voltage and frequency ratings of the power line are the same as those indicated on the unit's label. If this is not the case, do not connect the system to the power line until you adjust the unit to match the power source. In U.S.A, if the installation of this equipment will use 240V rather than 120V, the source must be a center-tapped, 240V, single-phase circuit.

Restricted Sale

U.S.A federal law restricts this device to sale by or on the order of a physician.

Supervised Use

This equipment is intended for use under the direct supervision of a licensed health care practitioner.

Ventilation Requirements

Set up the device in a location which affords sufficient ventilation. The ventilation openings of the device must not be obstructed. The ambient conditions specified in the technical specifications must be ensured at all times.

·Put the monitor in a location where you can easily see the screen and access the operating controls.

·This product is protected against the effects of cardiac defibrillator discharges to ensure proper recovery, as required by test standards. (the screen may blank during a defibrillator discharge but recovers within second as required by test standards.)

Reference Literature

Medical Device Directive 93/42/EEC

EN 60601-1/1990 +A1: 1993 +A2 : 1995 : Medical electrical equipment.

General requirements for safety

EN 60601-1-1/9. 1994 +A1 12.95: General requirements for safety.

General Precaution on Electric Safety

Warning

Check the items listed below before operating the equipment.

1. Be sure that AC power supply line is appropriate to use. (AC100 - 240V)
2. Be sure that the power adapter is the one supplied from Bionet. (DC18V, 2.5A)
3. Be sure that all cables are properly and firmly fixed.
4. Be sure that the equipment is properly grounded. (If not, this might cause a malfunction to occur in the product.)
5. The equipment should not be placed in the vicinity of electric generators, X-rays, broadcasting apparatus to eliminate the risk of electric noise during operation. This may cause incorrect results.

Note

The Equipment should be placed far from generator, X-ray equipment, broadcasting equipment or transmitting wires, so as to prevent the electrical noises from being generated during the operation, When these devices are near the Equipment, it can produce inaccurate measurements. For BM3Vet, both independent circuit and stable grounding are essentially required. In the event that the same power source is shared with other electronic equipment, it can also produce inaccurate output.

Warning

Do not operate the machine while in contact with the animal. It may cause serious danger to the user. Use only the provided patient cable.

Warning

In case the Equipment does not operate as usual or is damaged, do not use on animal, and immediately refer to the medical equipment technician of the hospital or the equipment supply division.

Note

BM3Vet is classified as follows:

- BM3Vet classifies as Class **I, BF & CF** concerning electric shock. It is not proper to operate this Equipment around combustible anesthetic or dissolvent.
- Noise level is B class regarding IEC/EN 60601-1 and the subject of Noise is B level concerning IEC/EN60601-1-2.

Equipment Connection

Caution

In the hospital, doctors and animals are exposed to dangerous, uncontrollable compensating currents. These currents are due to the potential differences between connected equipment. The safety solution to the problem is accomplished with EN60601-1;1993.

Biocompatibility

When used as intended, the parts of the product described in this operator manual, including accessories that come in contact with the animal during the intended use, fulfill the biocompatibility requirements of the applicable standards. If you have questions about this matter, please contact BIONET or its representatives.

Maintenance and Washing Equipment Connection

Various methods can be used to clean BM3Vet and its accessories. Please follow the methods mentioned below to avoid unnecessary damage or contamination to the Equipment.

We do not repair free of charge regardless of warranty period if it is contaminated or damaged due to use of dangerous materials not designated for washing.

Cleaning Applied Parts

Cables and Leadwires

CAUTION

Do not use acetone or ketone solvents for cleaning; do not use an autoclave or steam cleaner.

Cables and leadwires can be cleaned with a warm, damp cloth and mild soap, or isopropyl alcohol wipes. For more intensive disinfecting (near sterile) Ethylene Oxide (ETO) is acceptable but will reduce the useful lifetime of the cable or leadwire.

CAUTION

The decision to sterilize must be made per your institution's requirements with an awareness of the effect on the integrity of the cable or leadwire.

Note

The Equipment needs safety inspection once a year. Please refer to user's guide or service manual for the specifications.

Please check carefully both frame and sensor, after cleaning the Equipment, Do not use any equipment that is worn out or damaged.

At least once a month, clean and wipe off the frame by using the soft cloth after wetting it with water and alcohol. Do not use lacquer thinner, ethylene, or oxidizers which may lead to damage to the equipment.

Make sure both cables and accessories are free of dust or contaminants, and wipe them off with soft cloth wetted with warm water (40°), and at least once a week, clean them by using clinical alcohol.

Do not submerge the accessories under any liquid or detergent. Also, make sure any liquids do not

penetrate into the Equipment or probe.

Caution

Do not dispose of single use probe improperly. Always think about environmental contamination.

Caution

There is back-up battery on board inside system. When users dispose this battery, Please follow local laws for waste disposal.

Warning

Check the electrodes of batteries before changing them.

- Operate BM3Vet with internal electric power supply when unsure of external ground connection or installation occur.
- Remove the 1st Battery when not using equipment for an extended period of time.

For other applied parts such as temperature sensors, pulse oximetry probes, and NiBP cuffs, you must consult the manufacturer for cleaning, sterilization, or disinfecting methods.

1.3 Product Components

Product Outline

BM3Vet monitor is a product used for monitoring the biological vital signs of an animal. Main functions of the product include displaying information such as ECG, respiration, SpO2, NIBP and temperature on its LCD screen and monitoring parameter, and alarming. It also prints out waves and parameters via a printer.

Principal Characters of Product

BM3Vet is a compact, multifunctional monitoring equipment for an animal designed for portability. It features devices for auto power supply (DC 10V-16V) and DC power supply (DC 18V). The equipment also measures major parameters such as ECG, SpO2, NIBP, temperature and pulse, displaying it on a 7-inch color LCD screen. It also enables users to check waves and parameters and other vital signs of a animal via the 58mm thermal printer and monitor the animal by the remote-controlled alarm system. It also enables to build a central monitoring system by linking devices used for separate animals so that one can monitor several animals at a time.

Warning

Use only the supplement accessories provided by us. Otherwise, animal and user may be exposed to danger.

Warning

BEFORE USE — Before putting the system into operation visually inspect all connecting cables for signs of damage. Damaged cables and connectors must be replaced immediately. Before using the system, the operator must verify that it is in correct working order and operating condition. Periodically, and whenever the integrity of the product is in doubt, test all functions.

Product Configuration

1. Main body of BM3 Vet Monitor	1 EA
2. 3-Lead Animal Cable	1EA (3CBL-400, 3WIRE-400)
3. Disposable electrodes	10 EA (ECGSENS-400)
4. NIBP tubing (3M long)	1EA (NBPCBL-400)
5. Adult cuff (25-35 Cm)	1EA (ACUFF-400)
6. SpO2 sensor extension cable (2M)	1EA (SPCBL-400)
7. SpO2 Probe	1 EA (SPASENS-400)
8. DC Adaptor (MW160 made in AULT Co., Ltd.)	1 EA
9. Chart Paper (PAPER-400)	2ROLL

Option Product


1. Temperature (TEMPSENS-400)

Warning
In order to avoid electrical shock, do not open the cover. Disassembling of the equipment should be done only by the service personnel authorized by BIONET


Warning
Users must pay attention on connection any auxiliary device via LAN port or nurse calling. Always consider about summation of leakage current, please check if the auxiliary device is qualified by IEC 60601-1, or consult your hospital biomedical engineer.

Product Body Configuration


● FRONT

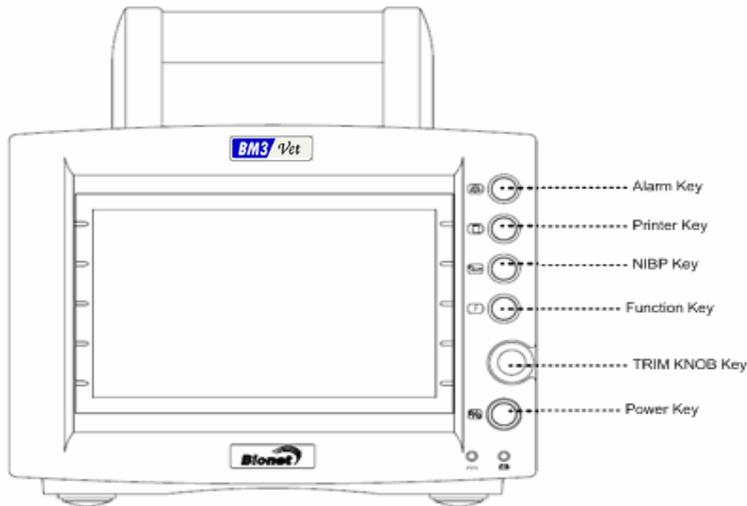


● BACK

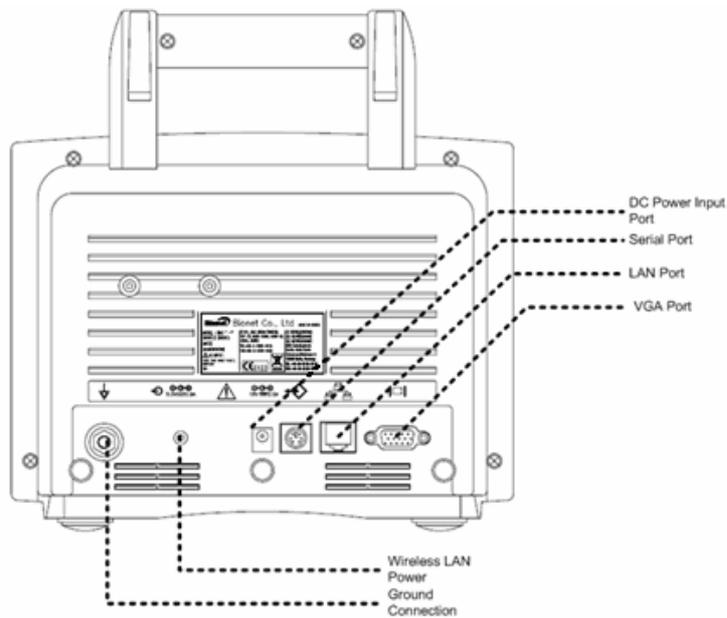

● Right Side

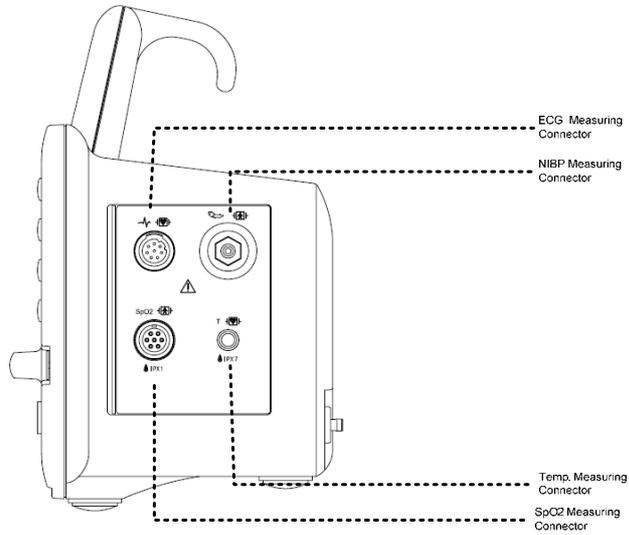

● Left Side

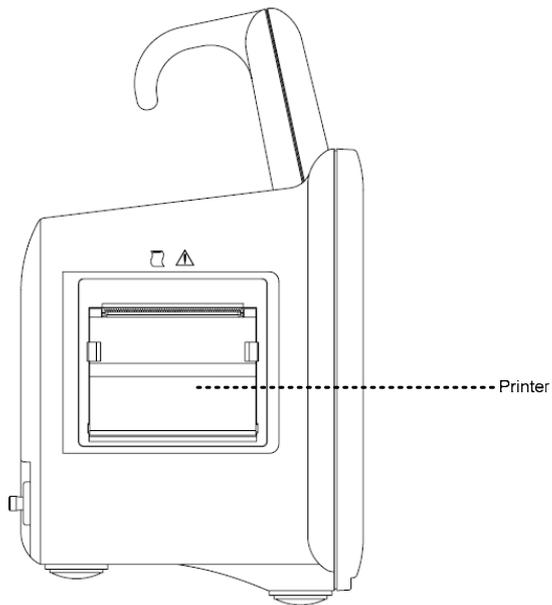


Accessories

ECG Cable +
Extension Cable



SpO₂ Cable +
Extension Cable

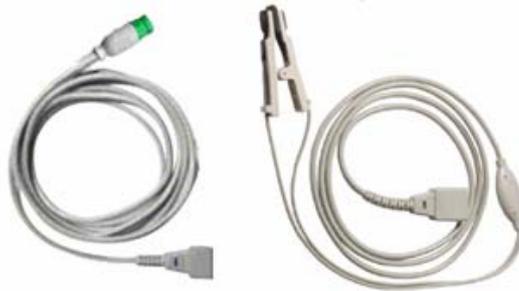

NIBP Cuff+
Extension cable






Temperature
sensor (Option)



Equipment Sign

	<p>ATTENTION :</p> <p>Consult accompanying documents</p>
	<p>TYPE CF APPLIED PART :</p> <p>Insulated (floating) applied part suitable for intentional external and internal application to the animal including direct cardiac application. "Paddles" outside the box indicate the applied part is defibrillator proof.</p> <p>Medical Standard Definition :</p> <p>F-type applied part(floating/insulated) complying with the specified requirements of IEC 60601-1/UL 2601-1/CSA 601.1</p> <p>Medical Standards to provide a higher degree of protection against electric shock than that provided by type CF applied parts.</p>
	<p>TYPE BF APPLIED PART :</p> <p>Insulated (floating) applied part suitable for intentional external and internal application to the animal excluding direct cardiac application. "Paddles" outside the box indicate the applied part is defibrillator proof.</p> <p>Medical Standard Definition :</p> <p>F-type applied part (floating/insulated) complying with the specified requirements of IEC 60601-1/UL 2601-1/CSA 601.1</p> <p>Medical Standards to provide a higher degree of protection against electric shock than that provided by type BF applied parts.</p>

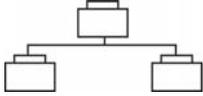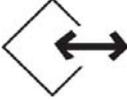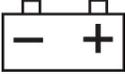	<p>External Ground</p>
	<p>Printer</p>
	<p>Serial Port</p>
	<p>LAN Port</p>
	<p>AUX Connector Port</p>
	<p>DC Input Indicator</p>
	<p>Battery Operation Indicator</p>
	<p>DC Input Connector</p>


	NIBP
T	Temperature
F	Function
	Power on
	Power off
	Respiration
	ECG
	Heart Pulse

1.4 Function and Key

External Function

The front panel of this product consists of an LCD screen and five function keys and one trim knob.


Operating the BM3Vet Patient Monitor

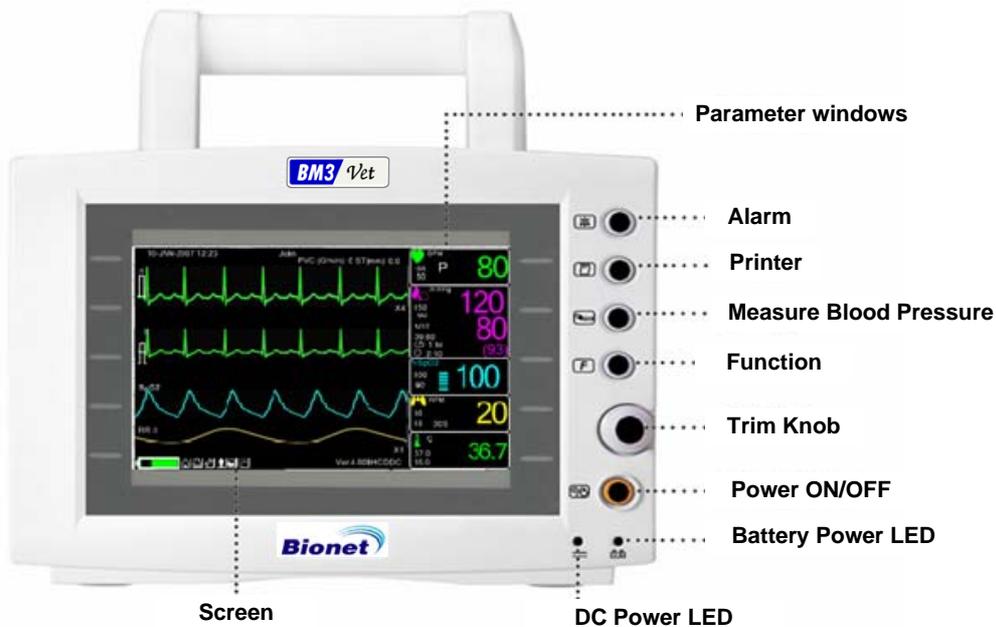


Operation Key

1. Power : Switches on and off the Power.
2. Function Key : Change Display mode.
3. Blood Pressure : Manually completes measuring blood pressure when this key is pressed with the cuff placed on the patient.
4. Printer : Prints out the waves selected from the menu until the key is pressed to stop.
5. Alarm : Stop alarm sound.
First press stops the current alarm for one minute
Second press stops the current alarm for five minutes.
Third press resets the alarm back to the original setting.
6. Trim Knob : This key is used to select menu options by turning it clockwise or counterclockwise to move cursors.

NIBP KEY

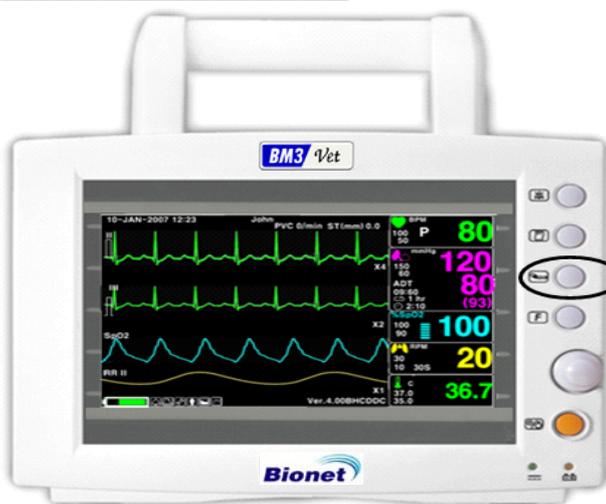

The manual press of external NIBP key can start or stop the operation of NIBP.

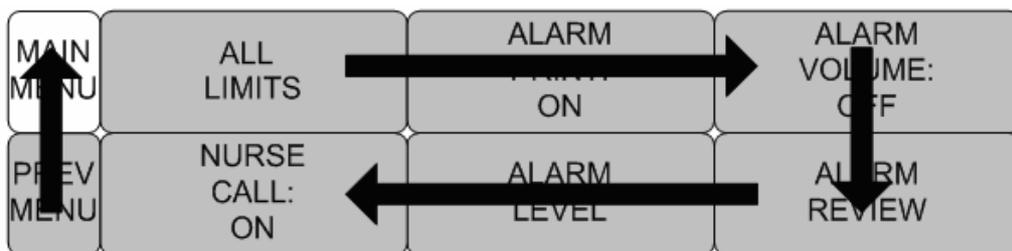


1.5 Standard Power Supply Application

DC Power

DC Power LED is lighted on when the DC Power is plugged into the inlet at the back of the product.

A press of the power key makes the machine ready for use.

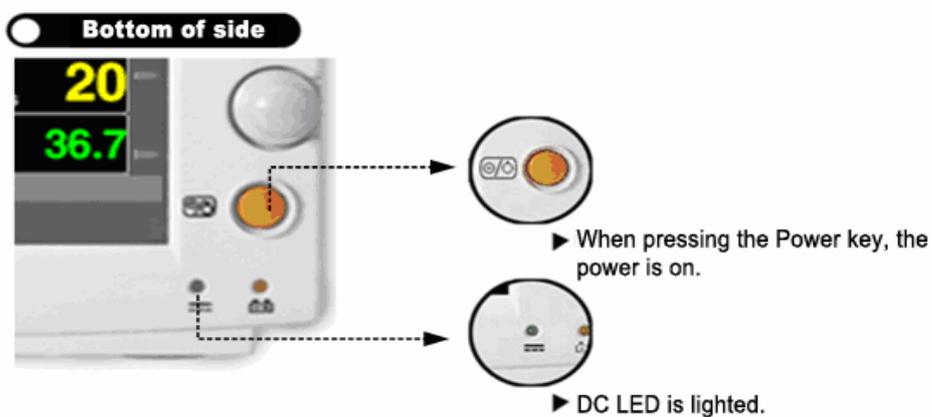

Warning

This equipment must only be connected to a grounded main.

1.6 Battery Power Supply Application

Battery power can be supplied for enabling portable use or use during AC power failure.

Operation

1. Battery Power LED is lit when the machine is in use.



2. The DC/battery power is only sustainable for 1 hour.

3. Battery is automatically charged when the machine is connected to DC Power Supply. Battery LED is lit on after blinking.

4. The charging status of the batteries is displayed with 5 green boxes, each indicating a different charging

. (0% -> 25% -> 50% -> 75% -> 100%)

- Battery: LS1865L2203S1PMXZ(11.1V - 2200mA, Li-ion)

The Lithium-Ion battery is a rechargeable battery containing Lithium-Ion cells. Each battery contains an integrated electronic fuel gauge and a safety protection circuit.



5. The discharge condition of battery is indicated with one of 5 yellow boxes, each box showing a different level of charge available.

. (100% -> 75% -> 50% -> 25% -> 0%)



When the battery power remains 25%, the message of “Low Battery” is displayed. The power is automatically cut off after 5 minutes from the appearance of the message. The machine will no longer operate when the “Low Battery” indication is on. Charge the batteries with the power adaptor, which BIONET provided.

-Battery charging time: More than 6 hours

-Continuous battery use time: Lowest 1 hour to highest 2 hours continuous use (buffering)

Warning
Check the electrodes of batteries before charging them.

6. Battery status indication: When battery is not installed or out of order, it is shown by a red 'X' as shown below.



7. Power from an automobile power adaptor: “CAR” appears rather than the battery symbols when an automobile power(12V~15V) is used.



Display of automobile power

Note
Battery is not charged when the automobile power is used.

The Impact of Lithium-Ion Battery Technology on the Battery

The following are the key points you should know about Lithium-Ion battery technology:

The battery will discharge on its own, even when it is not installed in a monitor. This discharge is the result of the Lithium-Ion cells and the bias current required for the integrated electronics.

By the nature of Lithium-Ion cells, the battery will self-discharge.

The self-discharge rate doubles for every 10°C (18°F) rise in temperature.

The capacity of the battery degrades significantly at higher temperatures.

As the battery ages, the full-charge capacity of the battery will degrade and be permanently lost. As a result, the amount of charge that is stored and available for use is reduced.

Conditioning Guideline

The battery in the monitor should be fully charged and discharged every six months.

Storage Guideline

Store the battery outside of the monitor at a temperature between 20°C to 25°C (68°F to 77°F).

When the battery is stored inside a monitor that is powered by an AC power source, the battery cell temperature increases by 15°C to 20°C (59°F to 68°F) above the room's ambient temperature. This reduces the life of the battery.

When the battery is stored inside a monitor that is continuously powered by an AC power source and is not powered by battery on a regular basis, the life of the battery may be less than 12 months. BIONET recommends that you remove the battery and store it near the monitor until it is needed for transport.


How to Recycle the Battery

When the battery no longer holds a charge, it should be replaced. The battery is recyclable. Remove the old battery from the monitor and follow your local recycling guidelines.

WARNING
EXPLOSION HAZARD —
DO NOT incinerate the battery or store at high temperatures. Serious injury or death could result.

1.7 General Manu Operation

Screen Composition

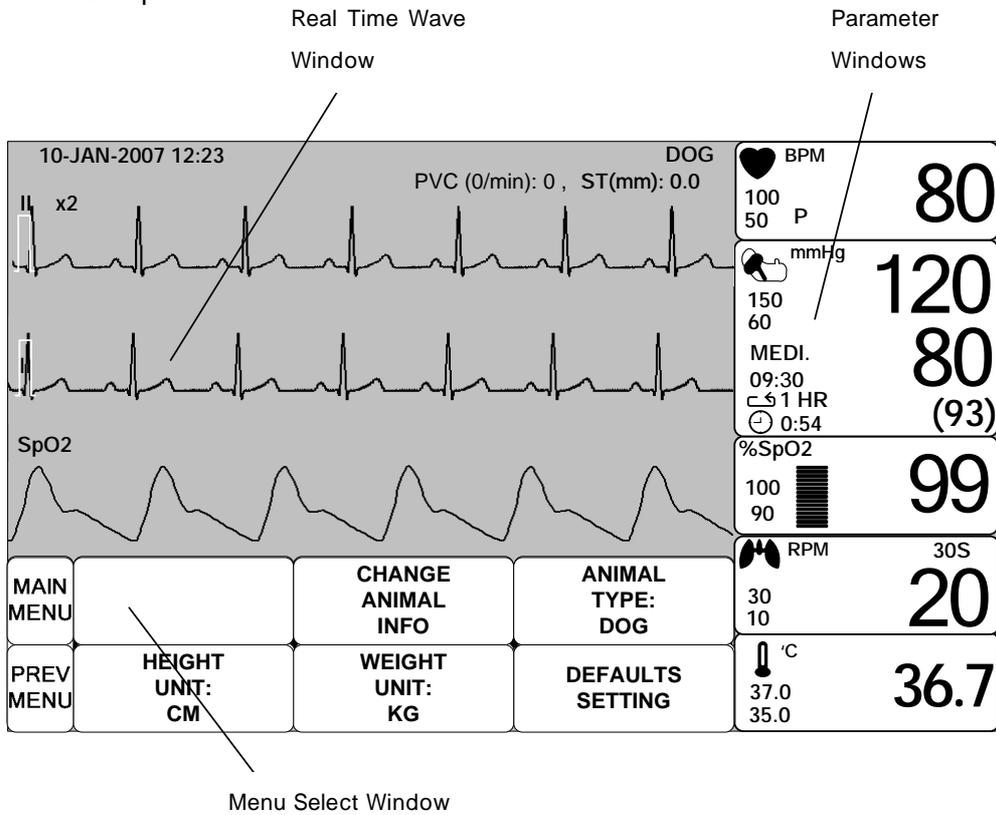


Real Time Wave Window : Displays measured results by up to three waves.

Menu Select Window : Menus appear when they are activated..

Parameter Window : Measured and setup data are displayed in five windows.

Menu Selection

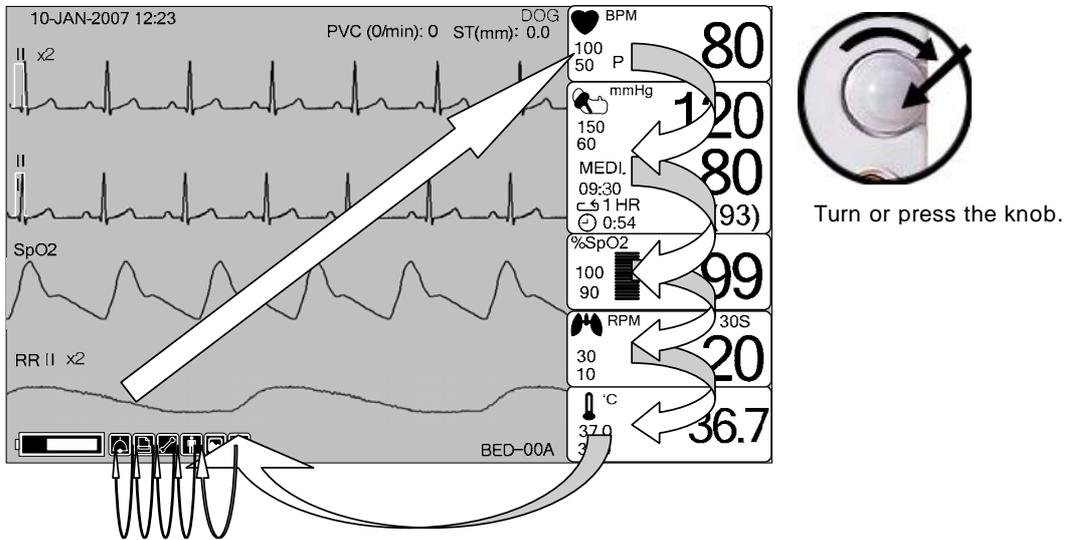

When the Trim Knob Key is turned, menus are selected in the order indicated above. The above screen shows that the MORE menu is selected. The menu move to the right in the order of MORE MENU → ECG → NIBP → SpO2 → RESP → TEMP.

Menu Composition


More Menu Window

When the additional menu is selected it will set and cancel the functions.

MAIN MENU	DISPLAY		USER SERVICE
PREV MENU	KEY SOUND: ON	DEMO: ON	MAKER SERVICE

Numerical value sign widow

This window displays a measured parameter, function setup, and the boundary of parameter values.

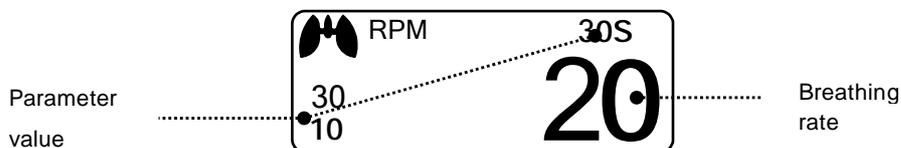

Menu selection by using Trim Knob key

As the key is turned to the right, the menu selection moves clockwise. As the key is turned to the left, the menu selection moves counterclockwise. The menu selection is activated when you depress Trim Knob key.

MAIN MENU	DISPLAY		USER SERVICE
PREV MENU	KEY SOUND: ON	DEMO : ON	MAKER SERVICE

Menu selection with arrows

Upward Movement: Turns the Trim Knob key to the left.

Downward Movement: Turns the Trim Knob key to the right.

Selection is made by pressing the Trim Knob key. One comes out of the menu after the selection.

MAIN MENU	LANGUAGE:	> ENGLISH FRENCH SPANISH ITALIAN GERMAN CHINESE
PREV MENU		

When moving within the quadrilateral, the letter reverses, and the numeric value reflects immediately.

MAIN MENU	QRS VOLUME : OFF	>	OFF	60%
			10%	70%
PREV MENU			20%	80%
			30%	90%
			40%	100%
			50%	

Word feature menu

The following figure shows the screen where the word sequence menu is activated within the word sequence correction menu. Here, the cursor moves over the words when the Trim Knob key is turned in the clockwise direction.

MAIN MENU	SET UNIT NAME	
PREV MENU		█

The above figure shows how the cursor moves on the screen. The cursor moves according to the direction the Trim Knob Key is turned. Press the Trim Knob key if you want to change a letter currently on the screen.


MAIN MENU	SET UNIT NAME	A
PREV MENU		█

The above figure shows how the cursor is selected to change a letter. Right-hand turning of the Trim Knob Key makes it possible to select in the order of 0-9,A-Z, and a blank, while left-hand turning makes the movement in the opposite direction. Once a letter or a number is selected, the screen comes back to the condition where the same process of selection can be made. One may move to

the menu item in the left of the screen to end the process, which is completed by pressing Trim Knob Key. After completion, the screen comes back to the earlier picture.

Operation menu

The setup value changes without a selection when the menu is moved.

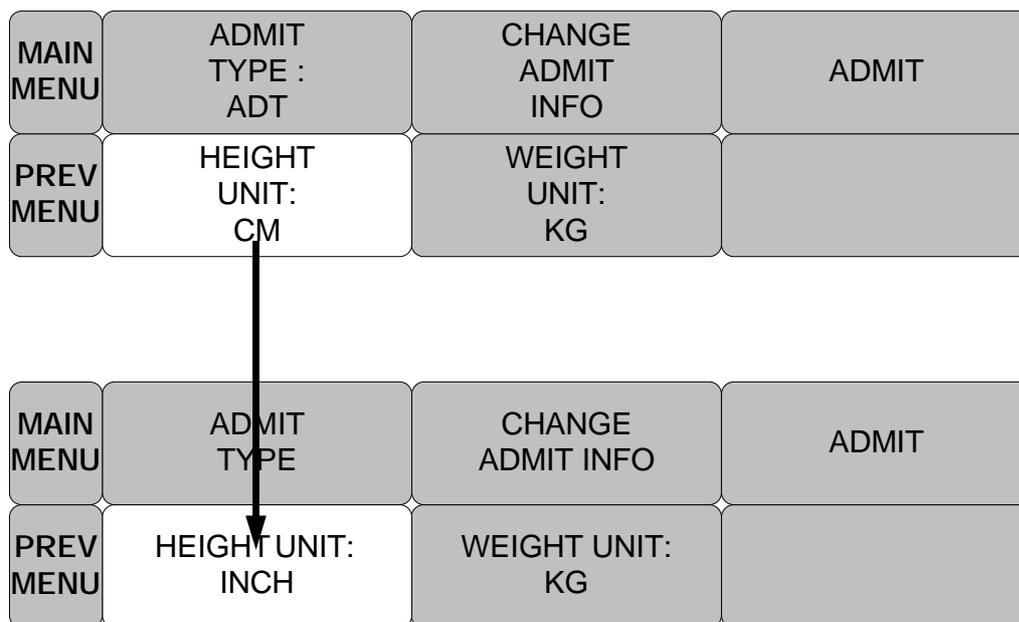

2. ANIMAL/DATA MANAGEMENT

2.1 ADMIT

CHANGE ANIMAL INFO
ANIMAL TYPE
HEIGHT
WEIGHT

2.2 ALARM

ALL LIMITS
ALARM PRINT
ALARM VOLUME
ALARM LEVEL
ARRHYTHM LEVEL
ALARM REVIEW
ALARM LIST
SAVE ALARM LEVEL
NURSE CALL

Additional setups are made for each parameter function. One can make an overall setup for the entire monitor system.

2.1 ADMIT

CHANGE ANIMAL INFO

ANIMAL TYPE

HEIGHT UNIT

WEIGHT UNIT

DEFAULTS SETTING

MAIN MENU		CHANGE ANIMAL INFO	ANIMAL TYPE: DOG
PREV MENU	HEIGHT UNIT: CM	WEIGHT UNIT: KG	

ANIMAL TYPE

You can select animal type as follow.

HORSE : LARGE ANIMAL // DOG: MEDIUM ANIMAL

PUPPY : SMALL ANIMAL // CAT : TINY ANIMAL

MAIN MENU		CHANGE ANIMAL INFO	ANIMAL TYPE: DOG
PREV MENU	HEIGHT UNIT: CM	WEIGHT UNIT: KG	DEFAULTS SETTING

MAIN MENU		ANIMAL TYPE: DOG	> HORSE DOG PUPPY CAT
PREV MENU	HEIGHT UNIT: CM		

CHANGE ANIMAL INFORMATION

Hospital ID(11 letters for each), animal name (11 letters for each), sex (male or female), date of birth, weight, height, and animal ID (11 characters)

MAIN MENU		CHANGE ANIMAL INFO	ANIMAL TYPE: DOG
PREV MENU	HEIGHT UNIT: CM	WEIGHT UNIT: KG	DEFAULTS SETTING

CHANGE ANIMAL INFORMATION	
> RETURN	CONTENTS
LAST NAME	
FIRST NAME	
ANIMAL ID	
SEX	MALE
BIRTH DATE	1 - JAN - 2000
AGE	0
HEIGHT	160.0 CM
WEIGHT	50.0KG

DEFAULTS SETTING

Animal information and all Alarm limits change to standard.

MAIN MENU		CHANGE ANIMAL INFO	ANIMAL TYPE: DOG
PREV MENU	HEIGHT UNIT: CM	WEIGHT UNIT: KG	DEFAULTS SETTING

HEIGHT

Unit of height is set as Cm / Inch.

MAIN MENU		CHANGE ANIMAL INFO	ANIMAL TYPE: DOG
PREV MENU	HEIGHT UNIT: CM	WEIGHT UNIT: KG	DEFAULTS SETTING

MAIN MENU		CHANGE ANIMAL INFO	ANIMAL TYPE: DOG
PREV MENU	HEIGHT UNIT: INCH	WEIGHT UNIT: KG	DEFAULTS SETTING

WEIGHT

Unit of weight is set as Kg / LBS.

MAIN MENU		CHANGE ANIMAL INFO	ANIMAL TYPE: DOG
PREV MENU	HEIGHT UNIT: CM	WEIGHT UNIT: KG	DEFAULTS SETTING

MAIN MENU		CHANGE ANIMAL INFO	ANIMAL TYPE: DOG
PREV MENU	HEIGHT UNIT: CM	WEIGHT UNIT: LBS	DEFAULTS SETTING

2.2 ALARM

Alarm is divided into two, alarm for the animal's condition and for the product's condition.

The animal's alarm sounds when the diagnostic functions (ASYSTOLE, VTAC/VFIB, and VTAC) are detected. Each alarm sound differs in order and volume according to the severity of HIGH, MEDIUM, LOW and MESSAGE.

HIGH		-5	≡	300	≡		
MEDIUM		-3	≡	300	≡		
LOW		-1	≡	300	≡		
MESSAGE			≡	300	≡		



: Alarm sounds



: Number flashes



: Waves are printed out



: Alarm lamp flashes

Alarm for the Product

The machine gives alarm sounds for its system with a related message flashing.

LOW  -1 ≡ Alarm Text ≡

ALARM LIMITS : The machine enables one to see and change the limits of alarm for all parameter functions.

ALARM PRINT : with an ON/OFF setup, the related information is printed out whenever an alarm is given.

ALARM VOLUME : volume of each alarm can be adjusted in 10 step.

ALARM LEVEL : Priority of each parameter alarm can be set up.

ALARM REVIEW : Shows the priority order information for all alarms of each measurement.

NURSE CALL : Set the ON/OFF feature of the NURSE CALL.

MAIN MENU	ALL LIMITS	ALARM PRINT: ON	ALARM VOLUME: OFF
PREV MENU	NURSE CALL: ON	ALARM LEVEL	ALARM REVIEW

It is able to see all the alarm range and change of measurement function.

ALL LIMITS

MAIN MENU	ALL LIMITS	ALARM PRINT: ON	ALARM VOLUME: OFF
PREV MENU	NURSE CALL: ON	ALARM LEVEL	ALARM REVIEW

ALL LIMITS			
RETURN	UNITS	LOW	HIGH
HR	BPM	50	150
SPO2-%	%	90	100
SPO2-R	BPM	50	150
RESP	RPM	10	30
RESP-A	SEC	0	20
NIBP-S	mmHg	80	200
NIBP-M	mmHg	60	140
NIBP-D	mmHg	20	120
TEMP	°C	30.0	42.0
ST	mm	-10.0	10.0
PVC	/min	0	20

ALARM PRINT

Set ON/OFF functions automatically. When the alarm is activated the corresponding information is printed on heat sensitive paper. Alarm level must be set to MEDIUM Level or above to print. But, LEAD FAULT AND LOW BATTERY Alarm do not trigger alarm print function.

MAIN MENU	ALL LIMITS	ALARM PRINT: ON	ALARM VOLUME: OFF
PREV MENU	NURSE CALL: ON	ALARM LEVEL	ALARM REVIEW

ALARM VOLUME

Set the alarm volume to be set at 10 grades.

MAIN MENU	ALL LIMITS	ALARM PRINT: OFF	ALARM VOLUME: OFF
PREV MENU	NURSE CALL: ON	ALARM LEVEL	ALARM REVIEW

MAIN MENU	ALARM VOLUME: OFF	> OFF	60%
		10%	70%
		20%	80%
PREV MENU		30%	90%
		40%	100%
		50%	

ALARM LEVEL

Set the order of priority in each alarm.

MAIN MENU	ALL LIMITS	ALARM PRINT: ON	ALARM VOLUME: OFF
PREV MENU	NURSE CALL: ON	ALARM LEVEL	ALARM REVIEW


MAIN MENU	PARAMETER LEVEL	ARRHYTHM LEVEL	
PREV MENU			

PARAMETER LEVEL

PARAMETER ALARM LEVELS	
RETURN	ALARM LEVEL
HR	MEDIUM
SPO2-%	MEDIUM
SPO2-R	LOW
RESP	MESSAGE
RESP-A	MESSAGE
NIBP-S	MEDIUM
NIBP-M	MEDIUM
NIBP-D	MEDIUM
TEMP	MESSAGE
LEAD FAULT	MESSAGE
LOW BATTERY	MESSAGE

ARRHYTH LEVEL

One can set up priorities when he or she uses the alarm for the diagnostic function.

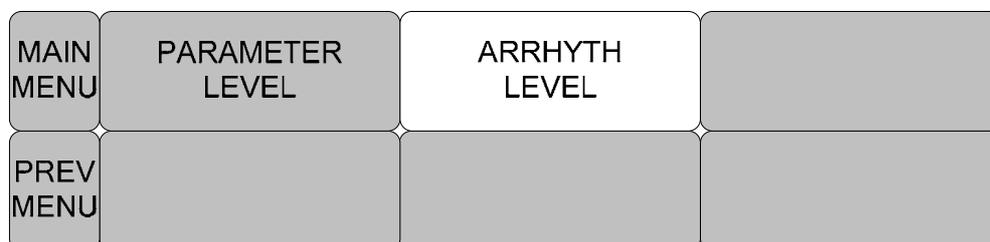

ARRHYTHMIA ALARM LEVELS	
RETURN	ALARM LEVEL
ASYSTOLE	HIGH
VTAC/VFIB	HIGH
VTAC	HIGH

ALARM REVIEW

After an alarm is triggered the alarms and data wave pattern can be reviewed. Set up for priority of each parameter alarm.

MAIN MENU	ALL LIMITS	ALARM PRINT: ON	ALARM VOLUME: OFF
PREV MENU	NURSE CALL: ON	ALARM LEVEL	ALARM REVIEW
MAIN MENU	ALARM LIST	SAVING CONDITION: HIGH	
PREV MENU			

ALARM LIST

When an alarm activates, this shows the order of the alarms.

MAIN MENU	ALARM LIST	SAVING CONDITION: HIGH	
PREV MENU			

10-MAR-2007 12:23
DOG

PVC (0/min): 0 ST(mm): 0.0

♥ BPM **80**
100
50 P

mmHg **120**
150
60

MEDI. **80**
09:30
1HR
0:53
(93)


%SpO2 **99**
100
90

👂 RPM **20**
30S
30
10

🌡 °C **36.7**
39.0
35.0


II x2

ALARM REVIEW		
RETURN	TIME	KIND
ECG	2007/03/10 10:22:45	HIGH
SPO2	2007/03/08 12:25:34	LOW
RESP	2007/03/06 23:32:10	HIGH
ECG	2007/03/05 09:12:36	MEDIUM
ECG	2007/03/04 13:52:42	MESSAGE
ECG	2007/03/03 18:18:38	MESSAGE
SPO2	2007/03/02 20:12:36	MESSAGE
SPO2	2007/03/01 22:25:56	MEDIUM
RESP	2007/03/01 09:12:15	MESSAGE
RESP	2007/02/26 14:52:38	MESSAGE
NIBP	2007/02/24 09:12:36	LOW

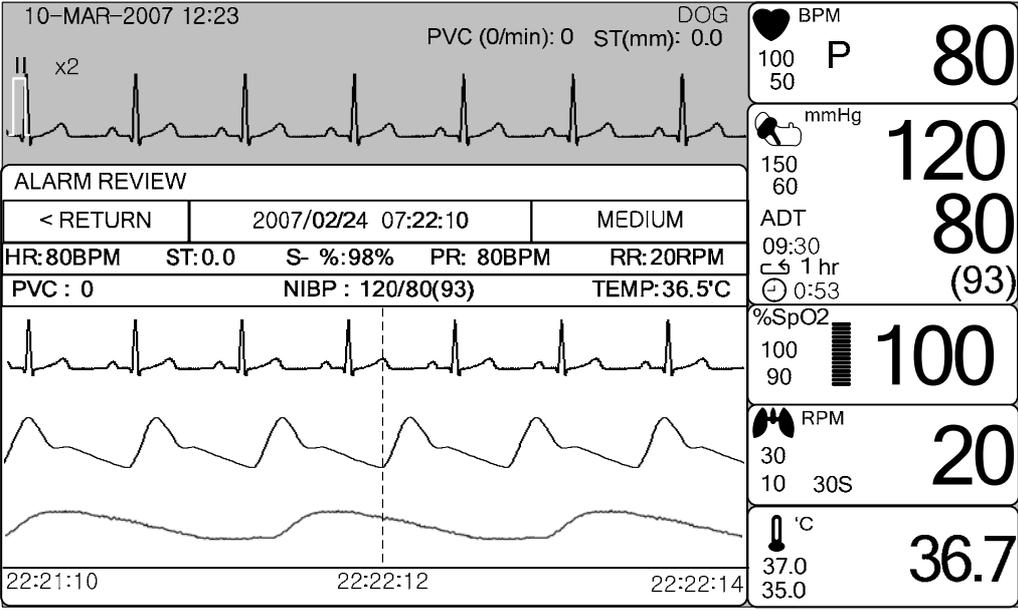

SAVING CONDITION

This determines the order in which triggered alarms are saved.

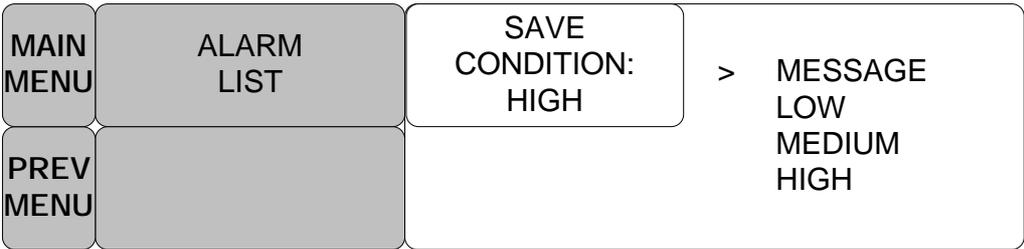

NURSE CALL

When an alarm is triggered, this activated the NURSE CALL function.

MAIN MENU	ALL LIMITS	ALARM PRINT: ON	ALARM VOLUME: OFF
PREV MENU	NURSE CALL: ON	ALARM LEVEL	ALARM REVIEW

MAIN MENU	ALL LIMITS	ALARM PRINT: ON	ALARM VOLUME: OFF
PREV MENU	NURSE CALL: OFF	ALARM LEVEL	ALARM REVIEW

3. SETUP

3.1 SETUP

DISPLAY

DEMO

USER SERVICE

MAKER SERVICE

3.1 SETUP

DISPLAY : screen set menu

USER SERVICE : This is the menu to set the connection used to interface with an external computer

MAKER SERVICE : This is the basic adjustment menu used to adjust the features of this product.

MAIN MENU	DISPLAY		USER SERVICE
PREV MENU	KEY SOUND: ON	DEMO: ON	MAKER SERVICE

DISPLAY

SET PARA : Measurement function selected.

WAVE SELECT : Set wave pattern source at the bottom of the WINDOW with LARGE PARAMETER MODE.

SET DATE & TIME: Set and change date and time.

HR SOURCE : Set and select HR/PR source.

COLOR SELECT: Set screen display color.

SET SWEEP: Set speed of ECG, RESP WAVE DISPLAY

MAIN MENU	SET PARA	WAVE SELECT: ECG	SET DATE & TIME
PREV MENU	SET SWEEP: 25mm/s		HR/PR SELECT: ECG

Commento [s1]: Wrong picture

SET PARA

Select measurement function to use

MAIN MENU	SET PARA	WAVE SELECT: ECG	SET DATE & TIME
PREV MENU	SWEEP SPEED: 25mm/s		HR SOURCE: ECG

PARAMETER WINDOW SET	
RETURN	WINDOW ON/OFF
ECG	ON
SPO2	ON
RESP	OFF
NIBP	OFF
TEMP	ON

WAVE SELECT

Select waveform to display in large parameter display.

MAIN MENU	SET PARA	WAVE SELECT: ECG	SET DATE & TIME
PREV MENU	SWEEP SPEED: 25mm/s		HR SOURCE: ECG

MAIN MENU	SET PARA	WAVE SELECT: ECG	> ECG SPO2 RESP
PREV MENU	SWEEP SPEED: 25mm/s		

SET DATE & TIME

It has sub menu to set date and time.

MAIN MENU	SET PARA	WAVE SELECT: ECG	SET DATE & TIME
PREV MENU	SWEEP SPEED: 25mm/s		HR SOURCE: ECG

SET TIME

Set time of equipment.

MAIN MENU	SET TIME	SET DATE	
PREV MENU			

MAIN MENU	SET TIME:	10 : 58 : 01
PREV MENU		

SET DATE

Set date of equipment

MAIN MENU	SET TIME	SET DATE	
PREV MENU			

MAIN MENU	SET DATE:	06-DEC-2007	
PREV MENU			

HR SOURCE

This menu is used to set the source that detects heart and pulse rate.

The source can select among ECG and SPO2.

MAIN MENU	SET PARA	WAVE SELECT: ECG	SET DATE & TIME
PREV MENU	SWEEP SPEED: 25mm/s		HR SOURCE: ECG

MAIN MENU	SET PARA	HR SOURCE: ECG	> ECG SPO2
PREV MENU	SWEEP SPEED: 25mm/s		

SWEEP SPEED

Set speed of drawing wave signal pattern in this widow.

MAIN MENU	SET PARA	WAVE SELECT: ECG	SET DATE & TIME
PREV MENU	SWEEP SPEED: 25mm/s		HR SOURCE: ECG

MAIN MENU	SWEEP SPEED: 25mm/s	> 6.25 mm/s 12.5 mm/s 25 mm/s 50 mm/s	SET DATE & TIME
PREV MENU			HR SOURCE: ECG

KEY SOUND

Set ON/OFF Key sound of equipment.

MAIN MENU	DISPLAY		USER SERVICE
PREV MENU	KEY SOUND: ON	DEMO: ON	MAKER SERVICE

DEMO

Set ON/OFF DEMONSTRATION of equipment.

MAIN MENU	DISPLAY		USER SERVICE
PREV MENU	KEY SOUND: ON	DEMO: ON	MAKER SERVICE

USER SERVICE

The user is able to set the communication parameters, power supply filter, and bed number.

MAIN MENU	DISPLAY		USER SERVICE
PREV MENU	KEY SOUND: ON	DEMO: ON	MAKER SERVICE

MAIN MENU	SET UNIT NAME	SET BED NUMBER : 00A	
PREV MENU	SYSTEM	AC FILTER: 50HZ	

SET UNIT NAME

Set up for Equipment name.

MAIN MENU	SET UNIT NAME	SET BED NUMBER : 00A	
PREV MENU	SYSTEM	AC FILTER: 50HZ	

MAIN MENU	SET UNIT NAME:	█	
PREV MENU			

SET BED NUMBER

Set up for animal bed number.

Allowable settings are from 0 to 9 and A to Z.

MAIN MENU	SET UNIT NAME	SET BED NUMBER : 00A	
PREV MENU	SYSTEM	AC FILTER: 50HZ	

MAIN MENU	SET UNIT NAME	SET BED NUMBER : 00A	0 0 A
PREV MENU	SYSTEM		

AC FILTER

AC FILTER is function where you can set power supply frequency. This feature is required because power supply frequency can be different from one country to another. . (The selectable frequencies are 50Hz and 60Hz.)

MAIN MENU	SET UNIT NAME	SET BED NUMBER : 00A	
PREV MENU	SYSTEM	AC FILTER: 50HZ	

MAIN MENU	SET UNIT NAME	SET BED NUMBER : 00A	
PREV MENU	SYSTEM	AC FILTER: 60HZ	

SYSTEM

System able to change and verify Equipment version information and system information

SYSTEM INFO SET	
RETURN	CONTENTS
MAIN VER	1.01.BVCDDCA
CENTRAL	ON
HOST IP	192 . 168 . 030 . 077
DEVICE IP	192 . 168 . 030 . 100
SUBNET	255 . 255 . 255 . 000
GATEWAY	192 . 168 . 030 . 001
MAC ADDR	00 : 02 : BD : 80 : CB : 00

MAKER SERVICE

Maker service is a menu is used by manufacturers.

MAIN MENU	DISPLAY		USER SERVICE
PREV MENU	KEY SOUND: ON	DEMO : ON	MAKER SERVICE

4. TREND

4.1 TREND

GRAPHIC TREND

TABLE TREND

TREND WINDOW SETUP

4.1 TREND

TREND shows saved data graphically displayed with numeric values.

Real-time data recording duration is 1 minute. Amount of saving time is for this data will be saving for 128hours.



 : Move to main screen

 : Move within the tables

 : Move up to other analysis function

 : Move down to other analysis function


 : Time(HOURS) period set menu at Graphic Trend

 : Time(MIN) period set menu at Tabular Trend

GRAPHIC TREND


Wave Data can be stored and seen according to section.

MAIN MENU	GRAPHIC TREND	TABULAR TREND	TREND WINDOW SETUP
PREV MENU			

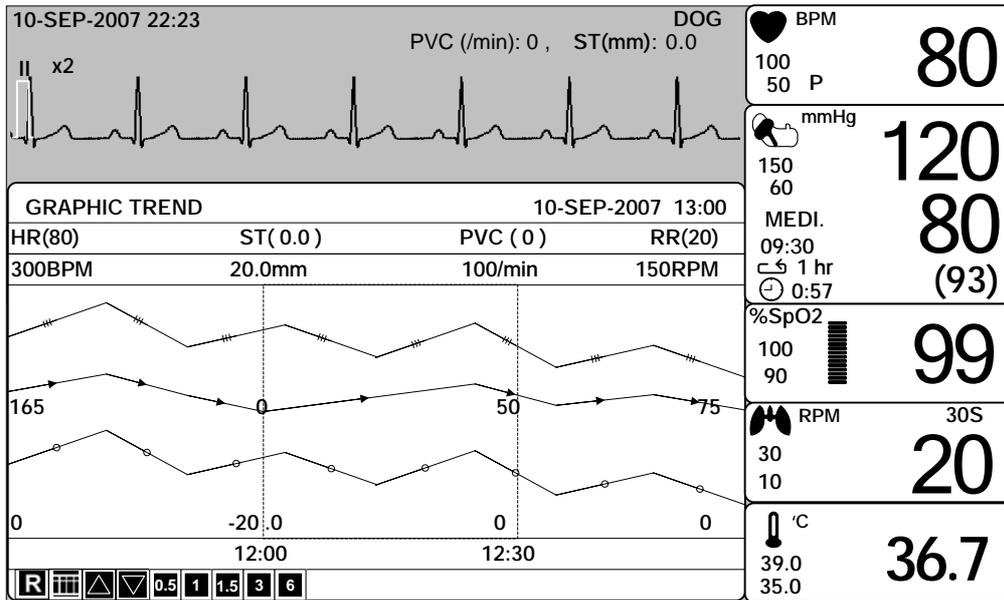

TIME PERIOD

One can set up and store data and time that one can see in a screen.

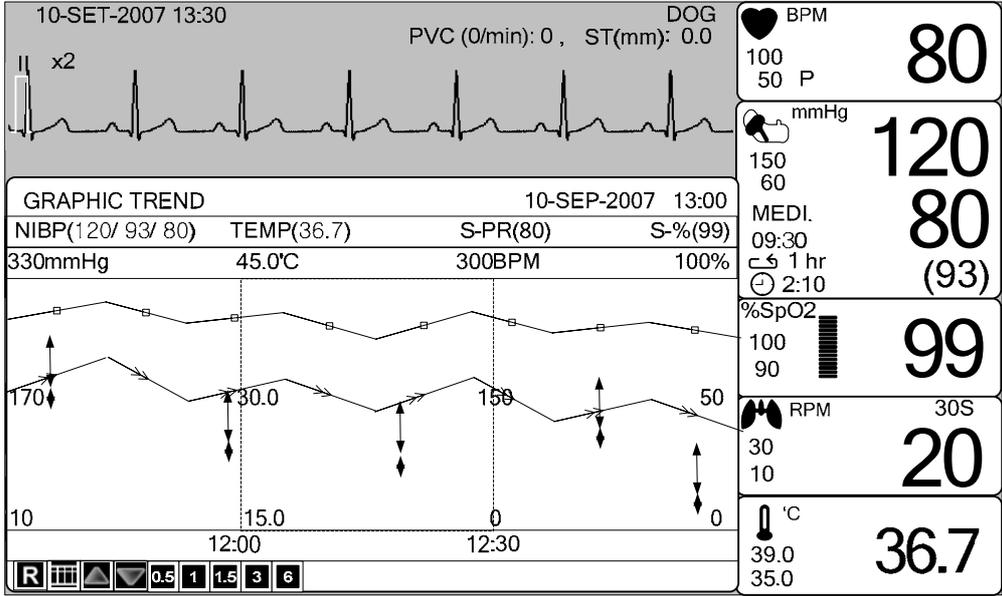

TABULAR TREND

One can see the stored data at the time previously set up.

MAIN MENU	GRAPHIC TREND	TABULAR TREND	TREND WINDOW SETUP
PREV MENU			

10-SEP-2007 13:30 DOG
PVC (0/min): 0 ST(mm): 0.0

BPM
100
50 P **80**

mmHg
150
60 **120**

MEDI.
09:30
1 hr
0:53 **80**
(93)

%SpO2
100
90 **99**

RPM
30S
30
10 **20**

°C
39.0
35.0 **36.7**

TABULAR TREND						10-SEP-2007 13:00	
	10-SEP 12:10	10-SEP 12:09	10-SEP 12:08	10-SEP 12:07	10-SEP 12:06		
HR	80	80	80	80	80		
SPO2-%	99	99	99	98	99		
SPO2-R	80	80	80	80	80		
RESP	20	20	20	20	20		
NIBP-S	120	120	120	120	120		
NIBP-M	93	93	93	93	93		
NIBP-D	80	80	80	80	80		
TEMP	36.7	36.7	36.7	36.7	36.7		
ST	0.0	0.0	0.0	0.0	0.0		
PVC	0	0	0	0	0		

R 1 5 15 30 60

/

TIME INTERVAL

One can store data and set up time.



10-SEP-2007 13:30 DOG

PVC (0/min): 0 ST(mm): 0.0

BPM 80

100 P
50

mmHg 120

150
60

MEDI. 80

09:30
1 hr
0:53

%SpO2 99

100
90

RPM 20

30S
30
10

°C 36.7


39.0
35.0

TABULAR TREND					
	10-SEP-2007 13:00				
	10-SEP 12:10	10-SEP 12:09	10-SEP 12:08	10-SEP 12:07	10-SEP 12:06
HR	80	80	80	80	80
SPO2-%	99	99	99	98	99
SPO2-R	80	80	80	80	80
RESP	20	20	20	20	20
NIBP-S	120	120	120	120	120
NIBP-M	93	93	93	93	93
NIBP-D	80	80	80	80	80
TEMP	36.7	36.7	36.7	36.7	36.7
ST	0.0	0.0	0.0	0.0	0.0
PVC	0	0	0	0	0

TREND WINDOW SETUP

Set the trend display window that will show the real time wave window.

MAIN MENU	GRAPHIC TREND	TABULAR TREND	TREND WINDOW SETUP
PREV MENU			

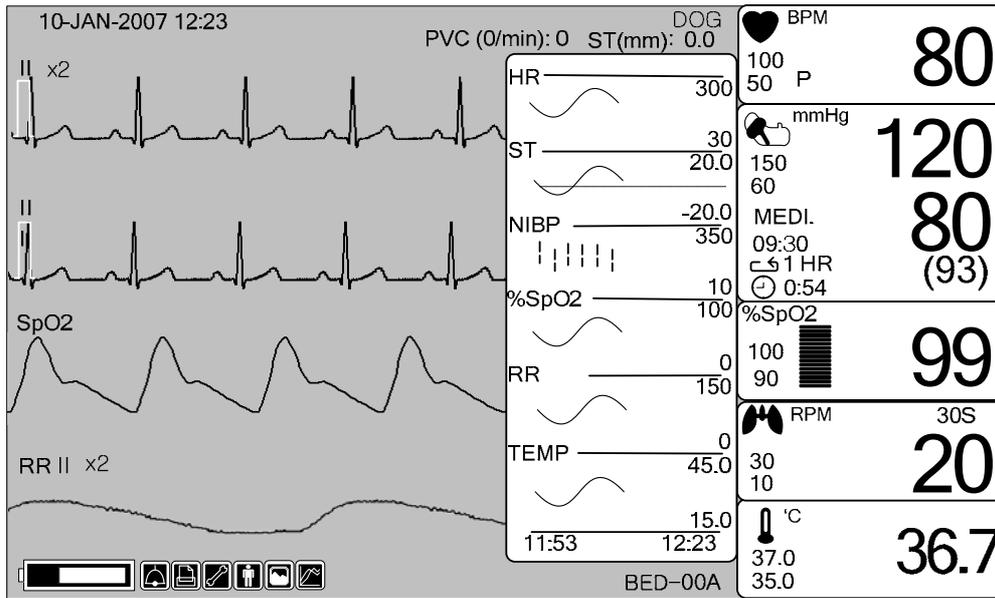

TIME PERIOD

Set visible time period in a screen.

MAIN MENU	TIME PERIOD: 30MINS	SET TREND PARA	
PREV MENU			

MAIN MENU	TIME PERIOD: 30MINS	> 30MINS 60MINS 90MINS 3HOUR 6HOUR 12HOUR	
PREV MENU			

SET TREND PARA

Set parameter for display in a screen.

MAIN MENU	TIME PERIOD: 30MINS	SET TREND PARA	
PREV MENU			

PARAMETER WINDOW SET	
RETURN	ON / OFF
HR	ON
ST	ON
SPO2	ON
PR	ON
RESP	ON
NIBP	ON
TEMP	ON

TREND PRINT

Graphic: After display a Graphic Trend on screen and press print to prints the selected trend.

Table: After display a Tabular Trend on screen and press print to receive print all the data in the selected animal admit (Admit) table.

5. ECG

5.1 Outline

Color and Name for Each Cable Size
ECG Connector Location and Measurement Cable
5 Lead Electrode Attached Location
3 Lead Electrode Attached Location
Method to Attach Electrode to Baby

5.2 ECG Data Window

5.3 ECG Data Setup

TRACE 1 LEAD SELECT
ALARM LIMIT
ALARM
QRS VOLUME
ECG SIZE
HEART RATE SOURCE
ECG SPEED
ANALYSIS SETTING

5.1 Introduction

It calculates the heart rate with 3 or 5 leads ECG signal acquisition and perform the alarm according to the setting value.

Colors and Standards of Cables


Leadwire	AHA Color code	AHA Label	IEC Color code	IEC Label
Right arm	White	RA	Red	R
Left arm	Black	LA	Yellow	L
Right leg	Green	RL	Black	N
Left leg	Red	LL	Green	F
V1(precordial)	Brown	V1	White	C1

AHA : American Heart Association (U.S.A. standard)

IEC : International Electro technical Commission (Europe standard)

Position of ECG Connector and Measuring Cable

ECG connector +detect cable

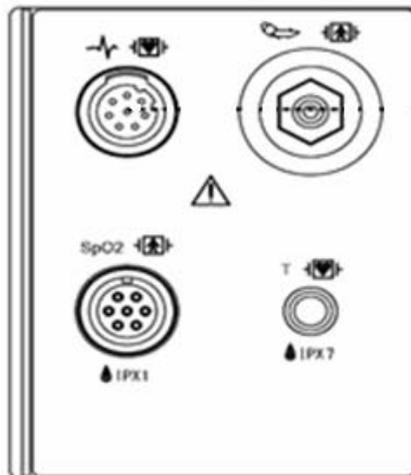



IEC 3LEAD CABLE



AHA 3LEAD CABLE



IEC 5LEAD CABLE



AHA 5LEAD CABLE



IEC 3LEAD



AHA 3LEAD



IEC 5LEAD



AHA 5LEAD

Attaching Electrodes to the Animal

1. Shave excess hair. With a piece of cotton pad moistened with alcohol, clean the animal's skin where the electrodes should be mounted. Avoid wrinkled or uneven skin areas. Wipe off the alcohol with a dry cotton pad.
2. Open the electrode package and take out the electrode.
3. Remove the backing paper from the electrode. Be careful not to touch the adhesive side.
4. Attach the disposable electrode to the previously cleaned skin. Avoid wrinkled and uneven skin areas.
5. The electrode lead which is connected to the monitor onto the electrode.
6. Fasten the electrode lead to the skin with surgical tape with an extra length of wire between the tape and the electrode. This prevents body movement from moving the electrode lead.

Note

- ✓ To maintain good contact between the electrode and skin, check that the paste of the disposable electrode is not dry.
- ✓ When contact of the disposable electrode becomes poor, replace the electrode with a new one immediately. Otherwise, contact impedance between the skin and electrode increase and the correct ECG cannot be obtained.
- ✓ If the contact is bad before the expiration date on the package, replace the electrode with a new one.
- ✓ To obtain a stable ECG waveform rub the skin with "skin Pure" skin preparation gel or tincture of Benzoin.
- ✓ Shall use only the CE certified disposable electrode.

Choosing an ECG lead for Arrhythmia Monitoring


It is very important to select a suitable lead for arrhythmia monitoring.

Guidelines for non-paced animals:

- ✓ QRS should be tall and narrow(recommended amplitude > 0.5mV)
- ✓ R wave should be above or below the baseline (but not bi-phasic)
- ✓ T wave should be smaller than 1/3 R-wave height.
- ✓ The P-wave should be smaller than 1/5 R-wave height.


To prevent detection of P-waves or baseline noises as QRS complexes, the minimum detection level for QRS complexes is set at 0.15mV. Adjusting the ECG wave size on the monitor display(gain adjustment)does not affect the ECG signal which is used for arrhythmia analysis. If the ECG signal is too small, you may get false alarms for asystole.

Information on the ECG waveform

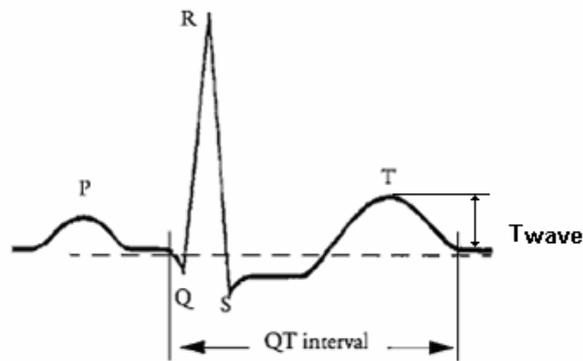


When ECG signal is 80bpm T-wave duration is 180ms, and the QT interval is 350ms.


Position of 3-Lead Wire Electrode

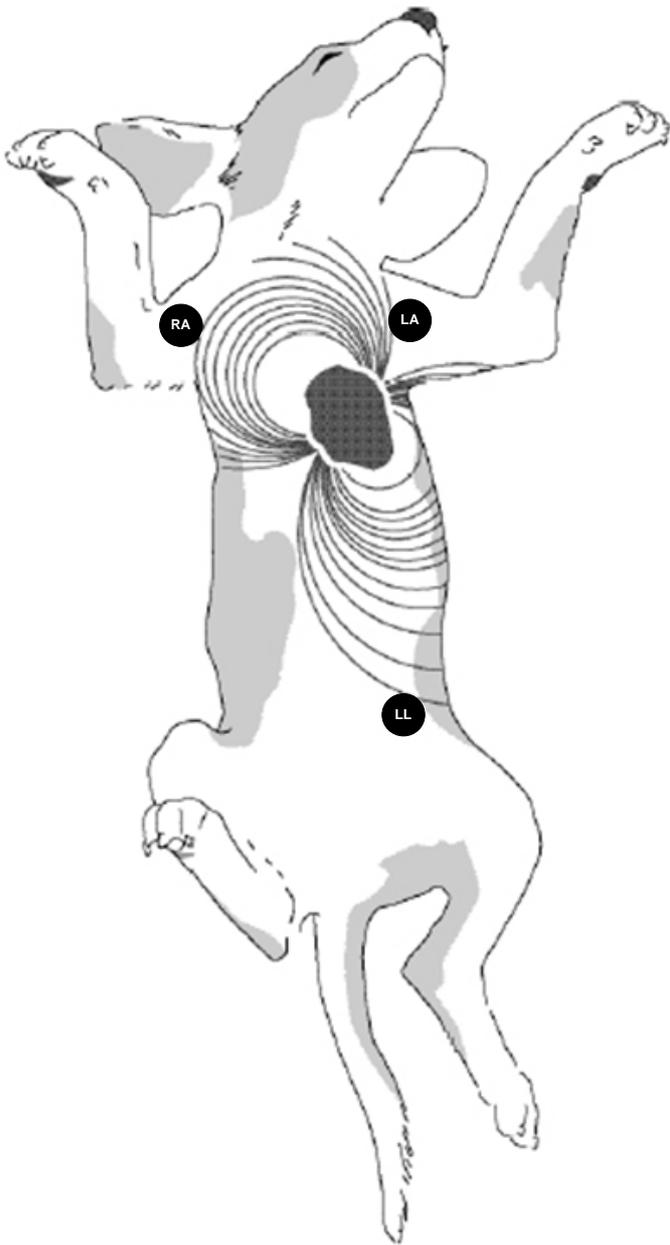

Position of 5-Lead Wire Electrode

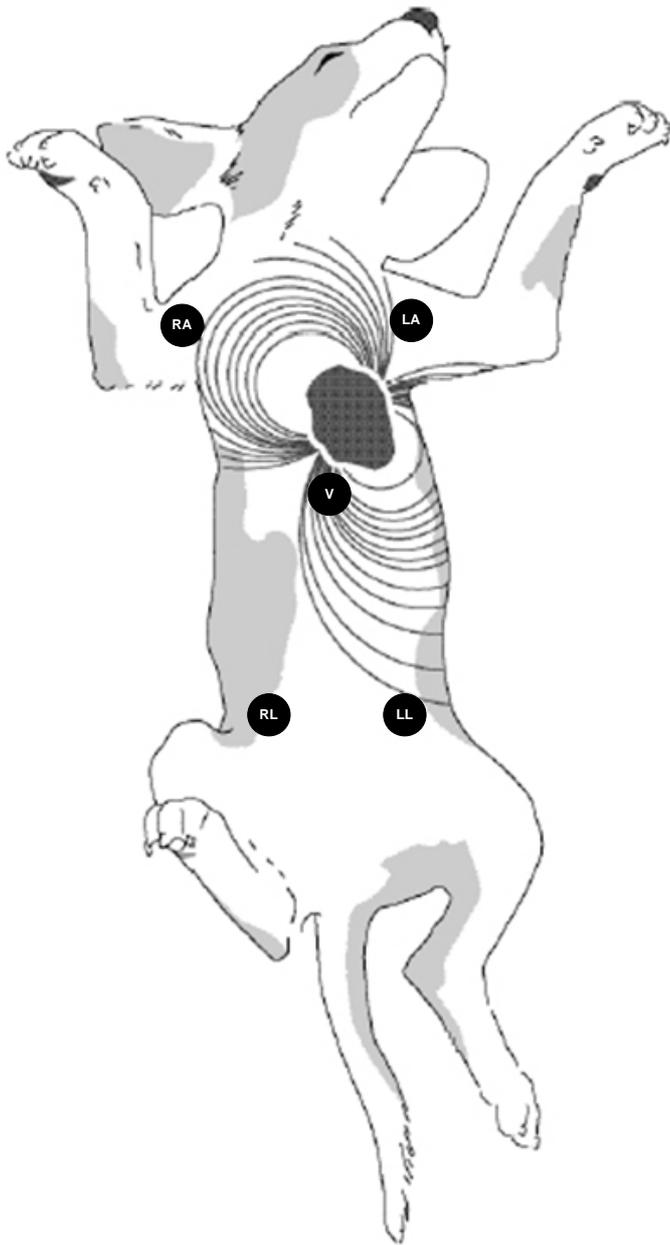

5.2 ECG Data Window

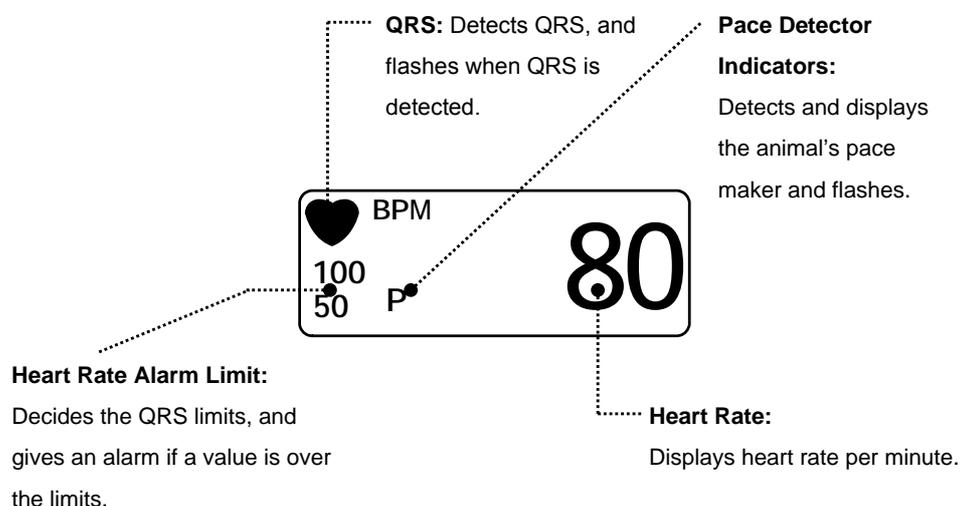

Note

ECG Wave Display is always on when the cable is connected.

The heart rate is calculated by a moving average. The monitor detects 8 consecutive beats, averages the R-R intervals of the latest 8 beats and uses this average to calculate the current heart rate. When a new beat is detected, the heart rate is recalculated using the latest 8 beats. The heart rate display is updated every 3 seconds.

Heart rate meter updates a new heart rate for a step increase or decrease in 10 seconds maximum.

When ventricular tachycardia is detected, the alarm set in 5 seconds maximum.

Check that the delay time of the output signal (alarm trigger 80ms maximum) is within the range of the connected equipment.

Safety Precautions

Warning

CABLES — Route all cables away from animal's throat to avoid possible strangulation.

CONDUCTIVE CONNECTIONS — Extreme care must be exercised when applying medical electrical equipment. Many parts of the animal/machine circuit are conductive, such as the animal, connectors, electrodes, transducers. It is very important that these conductive parts do not come into contact with other grounded, conductive parts when connected to the isolated animal input of the device. Such contact would bridge the animal's isolation and cancel the protection provided by the isolated input. In particular, there must be no contact of the neutral electrode and ground.

DEFIBRILLATION — Do not come into contact with animals during defibrillation. Otherwise serious injury or death could result.

To avoid the risk of serious electrical burn, shock, or other injury during defibrillation, all persons must keep clear of the bed and must not touch the animal or any equipment connected to the animal.

After defibrillation, the screen display recovers within 10seconds if the correct electrodes are used and applied in accordance with the manufacturer's instructions.

ECG cables can be damaged when connected to a animal during defibrillation. Check cables for functionality before using them again.

The peak of the synchronized defibrillator discharge should be delivered within 60ms of the peak of the R wave. The signal at the ECG output on the animal monitors is delayed by a maximum of 30ms.

If the ECG waveform on the screen is too unstable to synchronize with the animal's heart beat because of the following reason, remove the cause of an alarm, message, or unstable ECG, and then use a stable ECG lead for synchronization.

- ✓ ECG electrode is detached or broken. Lead wire is detached or broken.
- ✓ Lead wire moves. AC interference, EMG noise or noise from ESU is superimposed.
- ✓ Connection cable is broken or has a short circuit. Connector has poor contact.

INTERFACING OTHER EQUIPMENT — Devices may only be interconnected with each other or to parts of the system when it has been determined by qualified biomedical engineering personnel that there is no danger to the animal, the operator, or the environment as a result. In those instances where there is any element of doubt concerning the safety of connected devices, the user must contact the manufacturers concerned (or other informed experts) for proper use. In all cases, safe and proper operation should be verified with the applicable Manufacturer's instructions for use, and system standards IEC 60601-1-1/EN 60601-1-1 must be complied with.

Electrosurgery Unit

- ✓ Electrosurgical units(ESU) emit a lot of RF interference. If the monitor is used with an ESU,RF interference may affect the monitor operation.
- ✓ Locate the monitor as far as possible from the ESU. Locate them on opposite sides of the operating table, if possible.
- ✓ Connect the monitor and ESU to different AC outlets located as far as possible from each other.
- ✓ When using this monitor with an electrosurgical unit, its return plate and the electrodes for monitoring must be firmly attached to the animal. If the return plate is not attached correctly, it may burn the animal's skin where the electrodes are attached.

5.3 ECG Data Setup

A setup window appears at lower part of the screen when the Trim Knob Key is pressed in the ECG Parameter Window.

Selection is made by pressing the Trim Knob Key, while movement across the menu is performed by turning the key either clock or anticlockwise.

MAIN MENU	LEAD SELECT : II		ALARM
PREV MENU	DISPLAY	ANALYSIS SETTING	QRS VOLUME : OFF

TRACE 1 LEAD SELECT

Selection for display is made from channels I to V by moving the key left or right.

MAIN MENU	LEAD SELECT : II		ALARM
PREV MENU	DISPLAY	ANALYSIS SETTING	QRS VOLUME : OFF

MAIN MENU	LEAD SELECT : II	>	I	aVR
PREV MENU			II	aVL
			III	aVF
				V

ALARM LIMIT

Alarm Limit is 0 ~ 300.

1. Move the mark to select RETURN or HR, and press.
2. If pressed at HR, move to LOW, and press.
3. Checking the color changing, move right or left to select the value to set up, and press.
4. Press at HIGH. Checking the color changing, move again to select the value to set up, and press.
Move to HR, and press again. (You may decide to perform in the LOW HIGH order, which produces same result.)
5. Select RETURN to get out of the window.

MAIN MENU	LEAD SELECT : II		ALARM
PREV MENU	DISPLAY	ANALYSIS SETTING	QRS VOLUME : OFF
MAIN MENU	ALARM LIMIT	ALARM SOUND	
PREV MENU			

ECG ALARM LIMIT			
RETURN	UNITS	LOW	HIGH
HR	BPM	60	120

ALARM SOUND

Set ON/OFF of ECG alarm sound.

MAIN MENU	LEAD SELECT : II		ALARM
PREV MENU	DISPLAY	ANALYSIS SETTING	QRS VOLUME : OFF

MAIN MENU	ALARM LIMIT	ALARM SOUND	
PREV MENU			

ECG ALARM SOUND	
> RETURN	ECG ALARM SOUND
HR	ON
ARRHYTHMIA	ON
ST	ON
PVC	OFF

QRS VOLUME

Move the Key to select a volume rate from OFF, 10% to 100%.

MAIN MENU	LEAD SELECT : II		ALARM
PREV MENU	DISPLAY	ANALYSIS SETTING	QRS VOLUME : OFF

MAIN MENU	QRS VOLUME : OFF	>	OFF	60%
PREV MENU			10%	70%
			20%	80%
			30%	90%
			40%	100%
			50%	

DISPLAY

Set the sweep speed and waveform size.

MAIN MENU	LEAD SELECT : II		ALARM
PREV MENU	DISPLAY	ANALYSIS SETTING	QRS VOLUME : OFF

ECG SPEED

ECG speed is 25 mm/s.

Speed is changeable to 6.25, 12.5, 25, 50mm/s.

MAIN MENU	SWEEP SPEED : 25 mm/s	ECG SIZE : X1	HR SOURCE: ECG
PREV MENU			
MAIN MENU	SWEEP SPEED : 25 mm/s	6.25 mm/s 12.5 mm/s > 25 mm/s 50 mm/s	HR SOURCE: ECG
PREV MENU			

ECG SIZE

The size is changeable to X0.5, X1, X2, X4.

MAIN MENU	SWEEP SPEED : 25 mm/s	ECG SIZE : X1	HR SOURCE: ECG
PREV MENU			
MAIN MENU	SWEEP SPEED : 25 mm/s	ECG SIZE : X1	x 0.25 x 0.5 > x 1 x 2 x 4
PREV MENU			

HR SOURCE

MAIN MENU	SWEEP SPEED : 25 mm/s	ECG SIZE : X1	HR SOURCE: ECG
PREV MENU			

MAIN MENU	SWEEP SPEED : 25 mm/s	HR SOURCE: ECG	> ECG SPO2
PREV MENU			

ANALYSIS SETTING

Analysis setting divided to 3 menus.

ECG FILTER : One may select from three frequency types for WAVE FILTER.

MONITOR 0.5Hz ~ 40Hz

MODERATE 0.5Hz ~ 25Hz

MAXIMUM 5Hz ~ 25Hz

DIAGNOSIS 0.5Hz ~ 150Hz

MAIN MENU	LEAD SELECT : II		ALARM
PREV MENU	DISPLAY	ANALYSIS SETTING	QRS VOLUME : OFF

MAIN MENU	ECG FILTER : MONITOR	PACE : OFF	ARRHYTHM : OFF
PREV MENU		PVC SETTING	ST SETTING

MAIN MENU	ECG FILTER : MONITOR	> MONITOR MODERATE MAXIMUM DIAGONOSIS
PREV MENU		

PACE : Sets up ON/OFF to indicate that the animal has PACE.

The PACE menu option enables/disables the pacemaker detection program.

MAIN MENU	ECG FILTER : MONITOR	PACE : OFF	ARRHYTHM : OFF
PREV MENU		PVC SETTING	ST SETTING

Be aware of the following when monitoring a animal with a pacemaker.

ARRHYTHM : Sets up ON/OFF to indicate detection of diagnosis (ASYS, VTAC/VFIB, VTAC).

The Analysis algorithm simultaneously uses leads I, II, III, and the V lead for ECG and arrhythmia analysis.

MAIN MENU	ECG FILTER : MONITOR	PACE : OFF	ARRHYTHM : OFF
PREV MENU		PVC SETTING	ST SETTING

ASYSTOLE: Ventricular asystole occurs whenever the displayed heart rate drops to zero.

VTAC/VFIB: Ventricular fibrillation occurs when the ECG waveform indicates a chaotic ventricular rhythm with an average heart rate greater than or equal to 200beats per minute.

VTAC: Ventricular tachycardia occurs when a run of six or more ventricular beats is detected With an average heart rate greater than or equal to 150beats per minute.

ST SETTING : ST signal and settings related to ST menu.


MAIN MENU	ECG FILTER : MONITOR	PACE : OFF	ARRHYTHM : OFF
PREV MENU		PVC SETTING	ST SETTING

ST ANALYSIS: ON/OFF ST analysis signal.

MAIN MENU	ST ANALYSIS : ON	MEASUREMENT CONDITION	ST ALARM LIMIT
PREV MENU			ST ALARM LEVEL

MEASUREMENT CONDITION: ST measurement condition setting

MAIN MENU	ST ANALYSIS : ON	MEASUREMENT CONDITION	ST ALARM LIMIT
PREV MENU			ST ALARM LEVEL

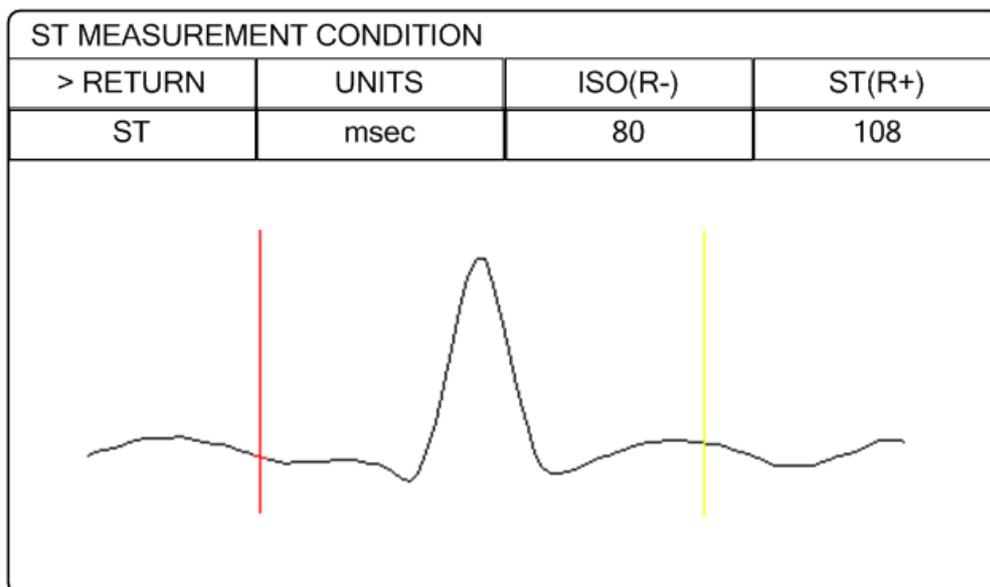

ST ALARM LIMIT: ST alarm limit range setting

MAIN MENU	ST ANALYSIS : ON	MEASUREMENT CONDITION	ST ALARM LIMIT
PREV MENU			ST ALARM LEVEL

ST ALARM LIMIT			
RETURN	UNITS	LOW	HIGH
ST	mm	-10.0	10.0

ST ALARM LEVEL: ALARM LEVEL setting

MAIN MENU	ST ANALYSIS : ON	MEASUREMENT CONDITION	ST ALARM LIMIT
PREV MENU			ST ALARM LEVEL

ST ALARM LEVEL	
> RETURN	ST ALARM LEVEL
ST	MEDIUM

PVC SETTING: PVC ON/OFF and ALARM limit range setting

MAIN MENU	ECG FILTER : MONITOR	PACE : OFF	ARRHYTHM : OFF
PREV MENU		PVC SETTING	ST SETTING

PVC ANALYSIS: Decision maker to display PVC value sign with ON/OFF

MAIN MENU	PVC ANALYSIS : ON		PVC ALARM LIMIT
PREV MENU			PVC ALARM LEVEL

PVC ALARM LIMIT: Set alarm indicate to PVC

MAIN MENU	PVC ANALYSIS : ON		PVC ALARM LIMIT
PREV MENU			PVC ALARM LEVEL

PVC ALARM LIMIT			
RETURN	UNITS	LOW	HIGH
PVC	/mm	0	20

PVC ALARM LEVEL: Set PVC ALARM LEVEL

MAIN MENU	PVC ANALYSIS : ON		PVC ALARM LIMIT
PREV MENU			PVC ALARM LEVEL

PVC ALARM LEVEL	
> RETURN	PVC ALARM LEVEL
PVC	MEDIUM

Warning
<p>Display Heart Beat Equipment Signal</p> <p>Heart Beat equipment signal displays when the PACE mode is on. The signal size or form are meaningless clinically.</p> <p>Number Of Heart Beats</p> <p>When using PACE mode, be sure to use an alternative method to monitor the heart rate. The heart beat equipment can show heart beat even during arrhythmia continuously. Therefore, do not depend on heart beat alarm excessively.</p>

CAUTION

FDA POSTMARKET SAFETY ALERT

The United States FDA Center for Device and Radiological Health issued a safety bulletin October 14, 1998. this bulletin states "that minute ventilation rate-adaptive implantable pacemakers can occasionally interact with certain cardiac monitoring and diagnostic programmed rate."

The FDA further recommends precautions to take into consideration for animals with these types of pacemakers. These precautions include disabling the rate responsive mode and enabling an alternate pace mode. For more information contact:

Office of Surveillance and Biometrics, CDRH, FDA
1350 Packard Drive, Mail Stop HFZ-510 Rockville, MD 20850 U.S.A

NOTE

ECG monitoring with animals in non-invasive trans coetaneous pacemakers may not be possible due to large amounts of energy produced by these devices. Monitoring ECG with an external device may be needed.

WARNINGS

VENTRICULAR ARRHYTHMISAS

The arrhythmia analysis program is intended to detect ventricular arrhythmias. It is not designed to detect a trial or supra ventricular arrhythmias. Occasionally it may incorrect identify the presence or absence of an arrhythmia. Therefore, a physician must analyze the arrhythmia information in conjunction with other clinical findings.

SUSPENDED ANALYSIS

Certain conditions suspend arrhythmia analysis. When suspended, arrhythmia conditions are not detected and alarms associated with arrhythmias do not occur. The messages which alert you to the conditions causing suspended arrhythmia analysis are : ARR OFF, ARRHSUSPEND, LEADS FAIL, ALARM PAUSE, ALL ALARMS OFF, and DISCHARGED.

FALSE CALLS

False low heart rate indicators or false asystole calls may result with certain pacemakers because of electrical overshoots.

MONITERING PACEMAKER PATIENTS

Monitoring of pacemaker patients can only occur with the pace program activated.

PACEMAKER SPIKE

An artificial pacemaker spike is displayed in place of the actual pacemaker spike. All pacemaker spikes appear uniform. Do not diagnostically interpret pacemaker spike size and shape.

PATIENT HAZARD

A pacemaker pulse can be counted as a QRS during a systole in either pace mode. Keep pacemaker patients under close observation.

RATE METERS

Keep pacemaker patients under close observation. Rate meters may continue to count the pacemaker rate during cardiac arrest and some arrhythmias. Therefore, do not rely entirely on rate meter alarms.

Trouble shooting

Problem :

Inaccurate heart rate and/or false asystole.

Solution :

Check ECG signal from animal:

1. Check/adjust lead placement.
2. Check/perform skin preparation.
3. Check/replace electrodes.

Check amplitude of ECG waveform:

1. Select ECG parameter label.
2. Select DISPLAY LEAD,
3. Scroll through all ECG leads and check for 0.5mV amplitude at normal (1X) size. (at least 0.5mV amplitude is required for QRS detection.) for borderline signals, validate on a graph.
4. If amplitudes are low, electrodes may need to be repositioned or replaced.

Problem :

False ventricular calls.

Solution :

Check ECG signal from animal: (the chest lead may exhibit polarity changes which may occasionally cause an inaccurate call.)

1. Check/adjust lead placement.
2. Check/perform skin preparation.
3. Check/replace electrodes. (if chest lead is a problem, move the chest lead to another chest position or leg position.)

Problem :

Inaccurate pacemaker detection

Solution :

Use pacemaker processing:

1. Select ECG parameter label.
2. Display the lead of ECG with the greatest amplitude in the top waveform position.
3. Select ANALYSIS SETTINGS.
4. SELECT DETECT PACE.

6. SpO₂

6.1 Outline

SpO₂ Connector Location and Measuring Cable

6.2 SpO₂ Data Window

6.3 SpO₂ Data Setup

SWEEP SPEED

RATE VOLUME

ALARM


ALARM LIMIT

6.1 Outline

SPO₂ monitoring is a noninvasive technique used to measure the amount of oxygenated hemoglobin and pulse rate by measuring the absorption of selected wavelengths of light. The light generated in the probe passes through the tissue and is converted into an electrical signal by the photodetector in the probe. The monitor processes the electrical signal and displays on the screen a waveform and digital values for SpO₂ and pulse rate. It detects SpO₂ in the way of transmitting the red and infrared rays into the capillary vessel to take the pulsation. Also performs the alarm function according to the setting values.

SpO₂ Connector Location and Measuring Cable


SpO₂ connector

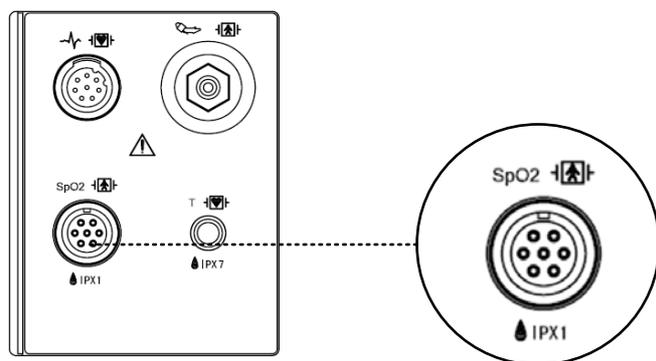

SpO₂ Cable +
Extension Cable



Position of SpO₂ Probe

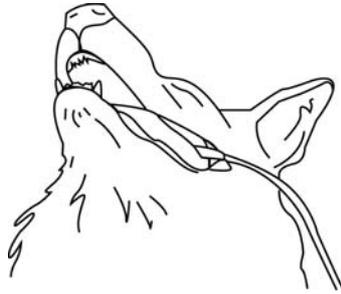


Note

The signal input is a high-insulation port and it is defibrillator proof ()

The insulated input ensures animal safety and protects the device during defibrillation and electro-surgery.

6.2 SpO₂ Data Window

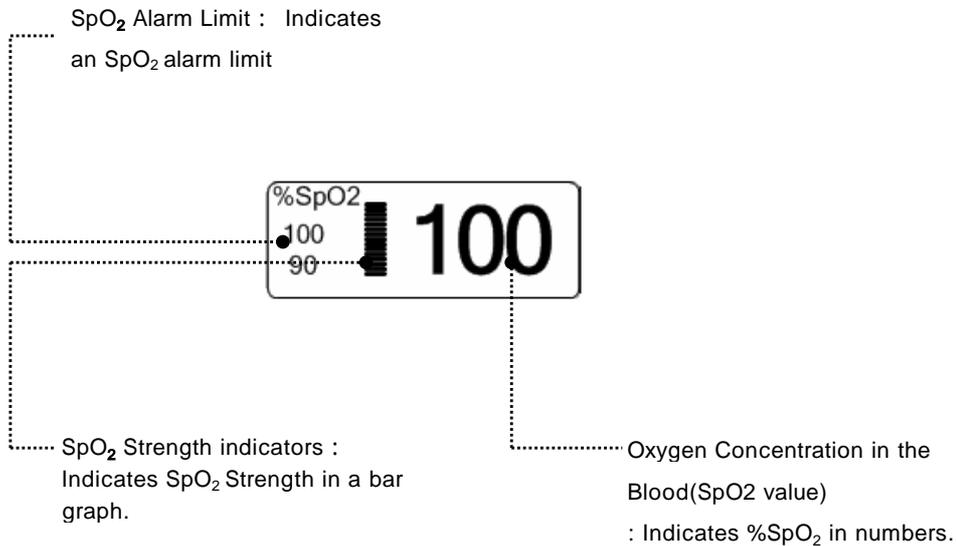

The current SPO2 value and the derived pulse rate (RATE) are displayed. The block sets indicate the strength of the signal (twenty block bars indicate the strongest signal). The SPO2 measurements are averaged over a 6-second period of time.

The monitor display is updated every second.

The SPO2 monitoring features are found in the SPO2 menu. These features include alarm limit adjustment, display of RATE, and RATE volume.

Note

SpO₂ WAVE SIZE is changed automatically.

Signal and Data Validity

It is extremely important to determine that the probe is attached to the animal correctly and the data is verifiable. To make this determination, three indications from the monitor are of assistance—signal strength bar, quality of the SPO2 waveform, and the stability of the SPO2 values. It is critical to observe all three indications simultaneously when ascertaining signal and data validity.

Signal Strength Bar

The signal strength bar is displayed within the SPO2 values window. This bar consists of 20 blocks set depending on the strength of the signal. Proper environmental conditions and probe attachment will help to ensure a strong signal.

Quality of SPO2 Waveform

Under normal conditions, the SPO2 waveform corresponds to (but is not proportional to) the arterial pressure waveform. The typical SPO2 waveform indicates not only a good waveform, but helps the user find a probe placement with the least noise spikes present. The figure below represents an SPO2 waveform of good quality.



Good Quality SPO2 Waveform

If noise (artifact) is seen on the waveform because of poor probe placement, the photodetector may not be flush with the tissue. Check that the probe is secured and the tissue sample is not too thick. Pulse rate is determined from the SPO2 waveform which can be disrupted by a cough or other hemodynamic pressure disturbances. Motion at the probe site is indicated by noise spikes in the normal waveform. (See the figure below.)



SPO2 Waveform with Artifact

Stability of SPO2 Values

The stability of the displayed SPO2 values can also be used as an indication of signal validity. Although stability is a relative term, with a small amount of practice one can get a good feeling for changes that are artifactual or physiological and the speed of each. Messages are provided in the SPO2 values window to aid you in successful SPO2 monitoring.

WARNING

In the monitoring of animals the coincidence of adverse conditions may lead to a disturbed signal going unnoticed. In this situation artifacts are capable of simulating a plausible parameter reading, so that the monitor fails to sound an alarm. In order to ensure reliable animal monitoring, the proper application of the probe and the signal quality must be checked at regular intervals.

6.3 SpO₂ Data Setup

ALARM : Menu in which SpO₂ Alarm are set up.

RATE VOLUME : Menu in which RATE VOLUME is set up

MAIN MENU	ALARM		RATE VOLUME: OFF

RATE VOLUME

Move the KEY to select the volume from OFF to 100%.

When the ECG volume rate is set, it turns OFF automatically.

MAIN MENU	ALARM		RATE VOLUME: OFF

MAIN MENU	RATE VOLUME: OFF	>	OFF	60%
			10%	70%
			20%	80%
			30%	90%
			40%	100%
			50%	

ALARM

Two menus: ALARM LIMIT, ALARM SOUND provided in the alarm menu

MAIN MENU	ALARM		RATE VOLUME: OFF

ALARM LIMIT

Number setting of alarm value of %SpO₂ is 0 ~ 100

1. Move the mark to select from RETURN, SpO₂ or SpO₂-R, and press.
2. After pressing at SpO₂, move the cursor right or left to LOW, and press.
3. Once the color is changed, move the cursor again to the selected value and press.
4. Place the cursor to HIGH and press, when the color changes, move the cursor again to select the targeted value, and press. Finally move to SpO₂ and press.
(You may decide to perform the process in the opposite order, LOW to HIGH, to have the same result.)
5. After pressing at SpO₂-R, move the cursor right or left to LOW, and press.
6. Once the color is changed, move the cursor again to the selected value and press.
7. Place the cursor to HIGH and press, when the color changes, move the cursor again to select the targeted value, and press. Finally move to SpO₂-R and press.
8. With the selection of RETURN the user gets out of the menu.

MAIN MENU	ALARM LIMIT	ALARM SOUND: ON	
PREV MENU			

SPO2 ALARM LIMIT			
RETURN	UNITS	LOW	HIGH
SPO2-%	%	90	100
SPO2-R	BPM	50	150

ALARM SOUND

Warning sound or message displays configuration menu when an alarm is triggered.

MAIN MENU	ALARM LIMIT	ALARM SOUND: ON	
PREV MENU			

MAIN MENU	ALARM LIMIT	ALARM SOUND: OFF	
PREV MENU			

LEAD FAULT Condition

When using a reusable lingual probe, there is a system alarm to alert you when the probe is off the Monitor. The monitor defaults this " LEAD FAULT" condition as a System Warning alarm. however, You can set it as a System ALARM LEVEL in Monitor Defaults.

SPO2 Messages

Below is a list of system status alarm messages which may be displayed in the SPO2 parameter window during monitoring.

CHECK PROBE

Reusable lingual probe is off the animal. Check the probe. *The factory default for this alarm is MESSAGE ALARM.*

PULSE SEARCH

Detection by the monitor of a repeatable pulse has ceased. Check the animal and the probe site.

POOR SIGNAL

The SPO2 signal is too low. No SPO2 data is displayed. This can be due to a low animal pulse, animal motion, or some other interference. Check the animal and the probe.

LOST SIGNAL

SPO2 data continues to be displayed, but the quality of the signal is questionable. Check the animal and the probe.

ARTIFACT

The SPO2 signal is patient's motion artifact and noise

7. RESPIRATION

7.1 Outline

Respiration Connector and Measuring Cable

7.2 RESPIRATION Data Window

7.3 RESPIRATION Data Setup

Respiration Size


Alarm Limit

7.1 Outline

Respiration via ECG Lead II electrode makes the skin area of the chest enlarged, causing changes in the resistance of skin. Through this it calculates respiration value per minute and performs the alarm function according to limit value.

Respiration Connector and Measuring Cable

Respiration Connector

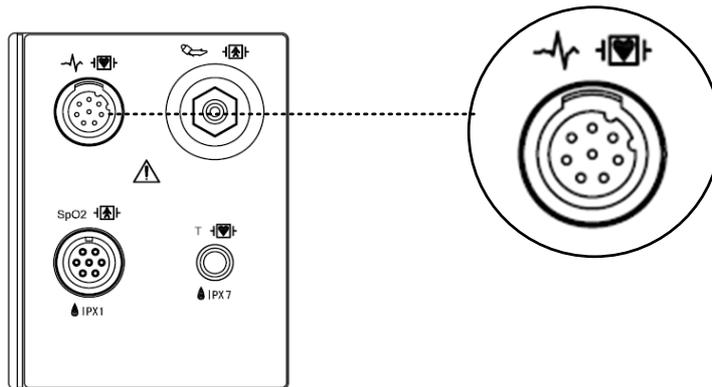

Respiration

Measuring Cable



IEC 3LEAD CABLE



AHA 3LEAD CABLE



IEC 5LEAD CABLE



AHA 5LEAD CABLE



IEC 3LEAD



AHA 3LEAD




IEC 5LEAD



AHA 5LEAD

7.2 Respiration Data Window

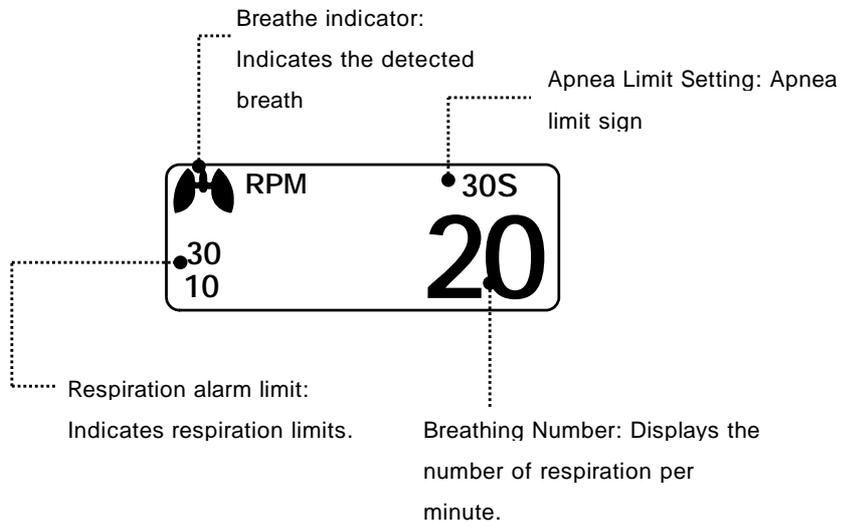

7.3 Respiration Data Setup

ALARM: Respiration alarm setting menu

RESP SIZE: A menu to setup Wave Display

SWEEP SPEED: A menu to setup Wave Display of speed

APNEA DETECT: A menu to setup APNEA alarm display

MAIN MENU	ALARM	SWEEP SPEED : 25mm/s	RESP SIZE : X 2
	APNEA DETECT : ON		

RESPIRATION SPEED

Wave pattern speed is 25 mm/s.

MAIN MENU	ALARM	SWEEP SPEED : 12.5mm/s	RESP SIZE : X 2
	APNEA DETECT : ON		

MAIN MENU	ALARM	SWEEP SPEED: 12.5mm/s	6.25 mm/s > 12.5 mm/s 25 mm/s
	APNEA DETECT : ON		

RESPIRATION

Set wave pattern size X2~ X10.

MAIN MENU	ALARM	SWEEP SPEED : 12.5mm/s	RESP SIZE : X 2
	APNEA DETECT : ON		

MAIN MENU	ALARM	RESP SIZE : X 2	> X 2 X 4 X 6 X 8 X10
	APNEA DETECT : ON		

APNEA DETECT

Deciding function of activating Apnea Alarm

MAIN MENU	ALARM	SWEEP SPEED : 12.5mm/s	RESP SIZE : X 2
	APNEA DETECT : ON		

MAIN MENU	ALARM	SWEEP SPEED : 12.5mm/s	RESP SIZE : X 2
	APNEA DETECT : OFF		

ALARM

Alarm menu provide ALARM LIMIT and ALARM.

MAIN MENU	ALARM	SWEEP SPEED : 12.5mm/s	RESP SIZE : X 2
	APNEA DETECT : ON		

ALARM LIMIT

Alarm Limit of Respiration Numeric Value is 5 ~ 150bpm

Alarm Limit of RESPIRATION APNEA Numeric Value is 3 ~ 30sec.

MAIN MENU	ALARM LIMIT	ALARM SOUND : ON	
PREV MENU			

1. Move the mark to select RETURN, RESP or RESP-A, and press.
2. After a press in RESP, move the cursor right or left to LOW, and press.
3. After the color changed, move the cursor right or left to the selected value, and press.
4. Place the cursor to HIGH, and press. When the color has changed, move the cursor again to select the value and press. Move to the RESP and press again. (You may decide to perform the process in the opposite order, LOW to HIGH, to have the same result.)
5. Once RESP-A is pressed, move to LOW and press.
6. When the color has changed, move the cursor to select the value, and press.
7. A press in the HIGH position, the color changes. Then move the cursor to select the value and press. Move again to RESP-A, and press.
8. Select RETURN to get out of the window.

RESP ALARM LIMIT			
RETURN	UNITS	LOW	HIGH
RESP	RPM	10	30
RESP-A	SEC	0	20

ALARM SOUND

Warning sound or message displays activation setting when Respiration ALRAM occurs.

MAIN MENU	ALARM LIMIT	ALARM SOUND : ON	
PREV MENU			

MAIN MENU	ALARM LIMIT	ALARM SOUND : OFF	
PREV MENU			

8. NIBP

8.1 Outline

NIBP Connector Location and Cuff

8.2 NIBP Data Window

8.3 NIBP Data Setup

ALARM LIMIT

ALARM

CUFF SIZE

UNIT SELECT

INTERVAL


INFLATION

8.1 Outline

This function is to measure minimum, Maximum and average blood pressure by using Oscillometric method

Position of NIBP Connector and cuff

NIBP Connector

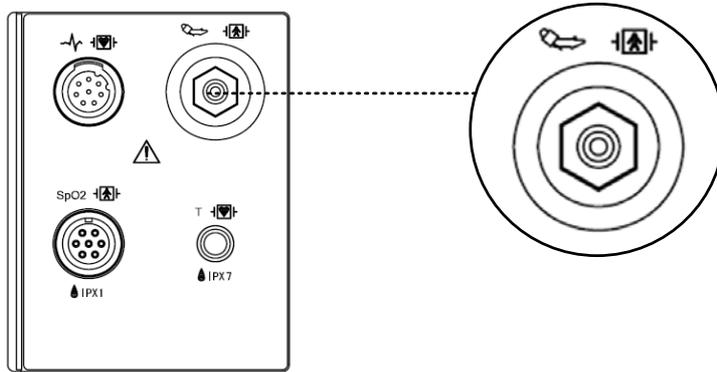

LARGE CUFF




Note

As the value of NIBP can vary according to the age and sex of an animal, the user needs to set up correct data in Parameter Menu before measurement.

WARNING

Noninvasive blood pressure monitoring is not recommended for animals with hypotension, hypertension, arrhythmias or extremely high or low heart rate. The software algorithm cannot accurately compute NIBP or animals with these conditions.

CAT CUFF Placement

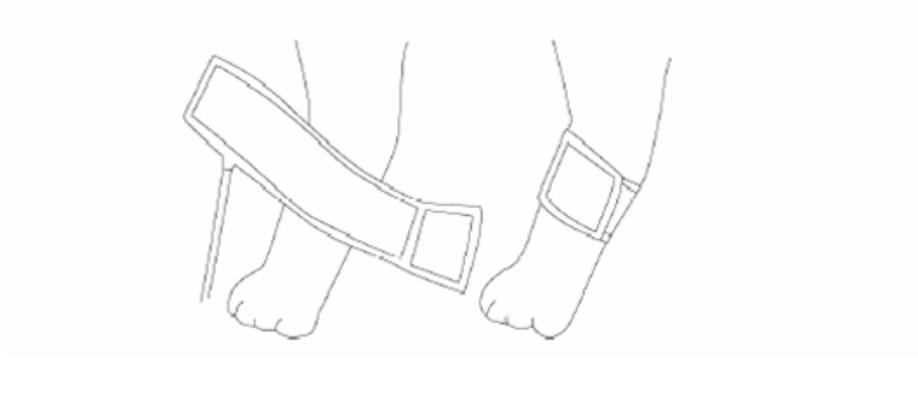

DOG CUFF Placement



Note

As the value of NIBP can vary according to the age and sex of a animal, the user needs to set up right data in parameter Menu before measurement. Tubes between the cuff and the monitor are not kinked or blocked.



The air pad should be exactly over the branchial artery. Tubing is immediately to the right or left of the branchial artery to prevent kinking when elbow is bent.

The routine maintenance is performed every 2 years.

Check the following list to ensure device operates properly and safety at all times.

1. Check for proper cuff size.
2. Check for residual air left in the cuff from a previous measurement.
3. Make sure cuff is not too tight or too loose.
4. Make sure cuff and heart are at same level, otherwise hydrostatic pressure will offset the NIBP value.
5. Minimize animal movement during measurement.
6. Check for leak in cuff or tubing.
7. Animal may have a weak pulse.

8.2 NIBP Data Window

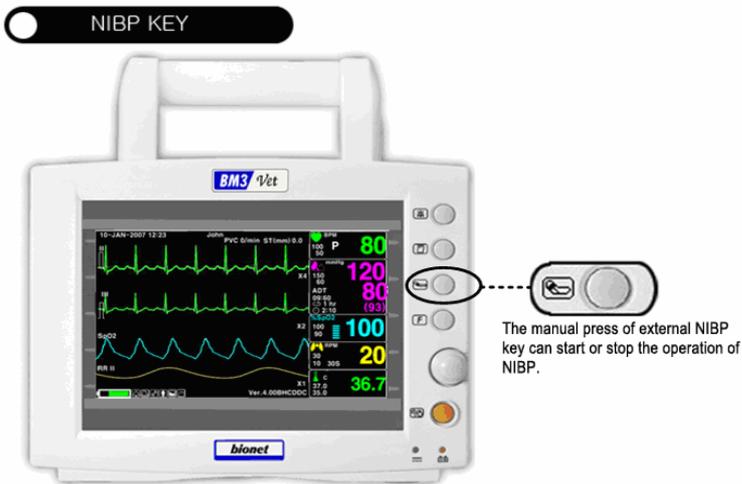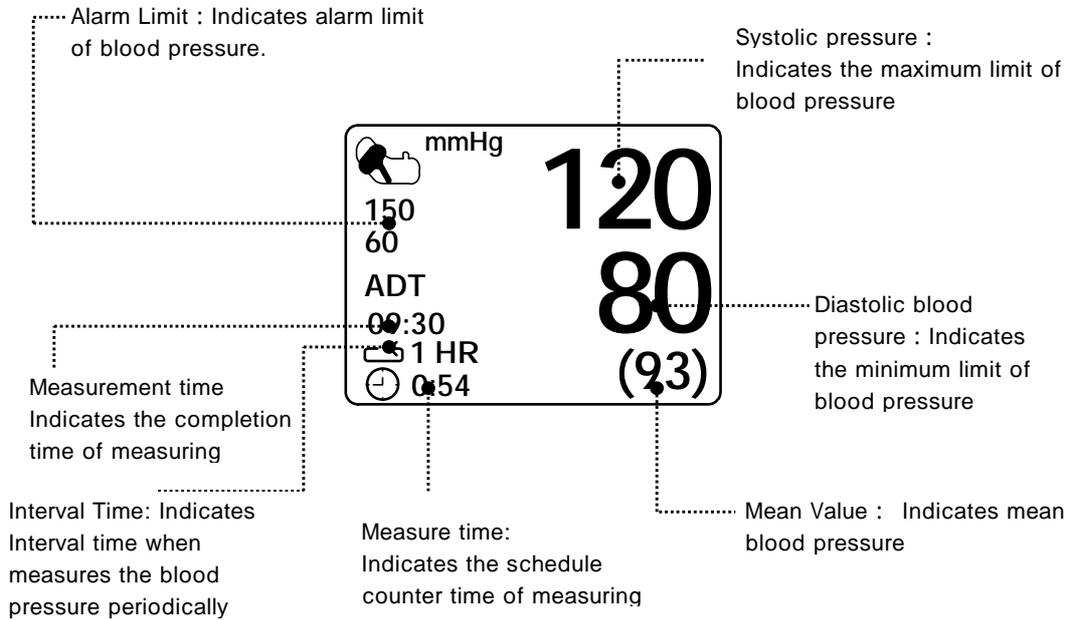

8.3 NIBP Data Setup

ALARM : A menu to set the Alarm

CUFF SIZE : A menu to select cuff size

UNIT SELECT: A menu to select the pressure unit

INTERVAL : A menu to set Interval time when measures the blood pressure periodically

INFLATION: Initial Pressurization setting menu

MAIN MENU	ALARM		CUFF SIZE: ADT
	UNIT SELECT: mmHg	INFLATION: 170mmHg	INTERVAL: OFF

ALARM

The alarm provides ALARM LIMIT and ALARM.

MAIN MENU	ALARM		CUFF SIZE: ADT
	UNIT SELECT: mmHg	INFLATION: 170mmHg	INTERVAL: OFF

MAIN MENU	ALARM LIMIT	ALARM SOUND: ON	
PREV MENU			

ALARM LIMIT

Alarm setting Numeric Value of Systolic, Diastolic, and mean pressure is 10 ~ 360mmHg.

1. Move the mark to select one from RETURN, NIBP-S, NIBP-M, or NIBP-D, and press.
2. Press the key at NIBP-S, and move to LOW, and press again. (The user gets the same result regardless of the LOW-HIGH, or HIGH-LOW order.)
3. When the color has changed, move it again to select a target value, and press.
4. Press the key at HIGH. When the color has changed, move to the right to select a target value, and press.
5. Set up or revise the values of NIBP-M and NIBP in the same way as above.
6. With the selection of RETURN, the user can get out of the window.

MAIN MENU	ALARM LIMIT	ALARM SOUND: ON	
PREV MENU			

NIBP ALARM LIMIT			
RETURN	UNITS	LOW	HIGH
NIBP-S	mmHg	80	200
NIBP-M	mmHg	40	140
NIBP-D	mmHg	20	120

ALARM SOUND

The menu which decides activation of warning signs and message displays when the respiration alarm is on.

MAIN MENU	ALARM LIMIT	ALARM SOUND: ON	
PREV MENU			

CUFF SIZE

The user can select a CUFF between LARGE and SMALL.

MAIN MENU	ALARM		CUFF SIZE: MEDIUM
PREV MENU	UNIT SELECT: mmHg	INFLATION: 170mmHg	INTERNAL: OFF

MAIN MENU	ALARM	CUFF SIZE: MEDIUM	> LARGE MEDIUM SMALL
PREV MENU	UNIT SELECT: mmHg		

UNIT SELECT

It is a function to set blood pressure measurement unit.

The blood pressure measurement unit provides mmHg and kPa.

MAIN MENU	ALARM		CUFF SIZE: ADT
	UNIT SELECT: mmHg	INFLATION: 170mmHg	INTERVAL: OFF

MAIN MENU	ALARM		CUFF SIZE: ADT
	UNIT SELECT: kPa	INFLATION: 170mmHg	INTERVAL: OFF

INTERVAL

This menu is used for selecting intervals when measures the blood pressure automatically.

Select a target interval from 1min, 2, 3, 4, 5, 10, 15, 20, 30, 1hour, 2, 4, 8.

MAIN MENU	ALARM		CUFF SIZE: ADT
	UNIT SELECT: mmHg	INFLATION: 170mmHg	INTERVAL: OFF
MAIN MENU	INTERVAL: OFF	> OFF 15MINS 1MIN 20MINS 2MINS 30MINS 3MINS 1HR 4MINS 2HRS 5MINS 4HRS 10MINS 8HRS	

INFLATION

It is a function for pressurization pressure.

Numeric value is 80 ~ 240 (10mmHg / step)

MAIN MENU	ALARM		CUFF SIZE: MEDIUM
PREV MENU	UNIT SELECT: kPa	INFLATION: 140mmHg	INTERVAL: OFF

MAIN MENU	ALARM		CUFF SIZE: MEDIUM
	UNIT SELECT: mmHg	INFLATION: 80mmHg	INTERVAL: OFF

MAIN MENU	ALARM		CUFF SIZE: MEDIUM
	UNIT SELECT: mmHg	INFLATION: 240mmHg	INTERVAL: OFF

Warning

Pay attention to not to block connecting hose when you put cuff on animal.

9. TEMPERATURE

9.1 Outline

Temperature Connector and Measuring Cable

9.2 Temperature Data Window

9.3 Temperature Data Setup

ALARM LIMIT


UNIT SELECT

9.1 Outline

This function is used to indicate the changes of resistance generated by the changes of temperature in numbers. The function involves the process of transferring the changes into electric signals.

Temperature Connector and Measuring Cable


Temperature Connector

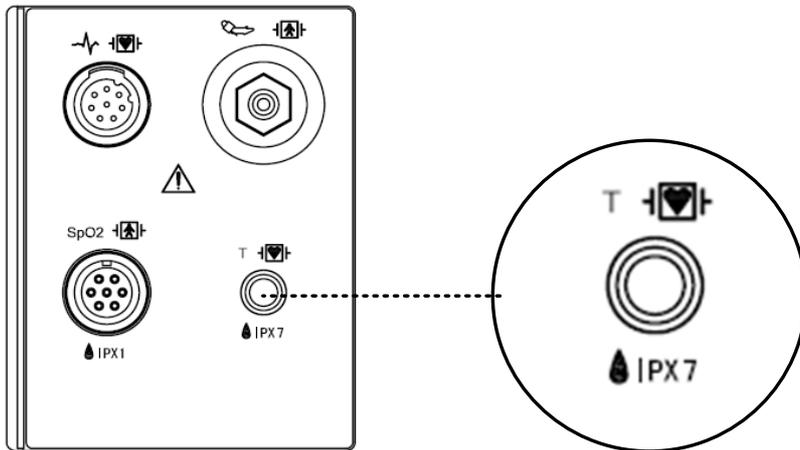

Temperature
Measuring Cable



9.2 Temperature Data Window

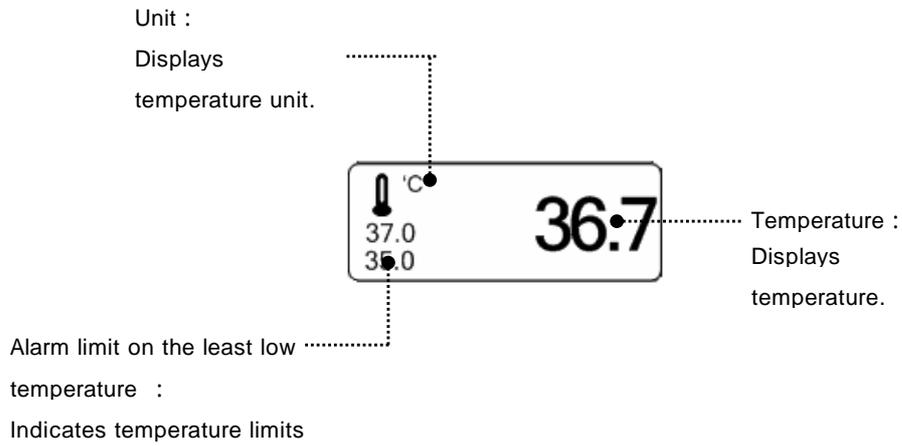

Note

The minimum measuring time required to obtain accurate readings at the specific animal site is at least 3 minutes.

9.3 Temperature Data Setup

ALARM : Temperature measurement alarm set

UNIT: Temperature measurement unit set

MAIN MENU	ALARM		UNIT SELECT: °C

ALARM

Alarm menu provide ALARM LIMIT and ALARM.

MAIN MENU	ALARM LIMIT	ALARM SOUND : ON	
PREV MENU			

ALARM LIMIT

Setting numeric value is 15.0 ~ 45.0 .

1. Move the mark to select either RETURN or TEMP, and press.
2. After pressing the cursor at TEMP, move it to LOW, and press.
3. When the color has changed, move the cursor again to select a target value, and press.
4. Move the cursor to HIGH and press. After the color has changed, move the cursor again to select a target value, and press. (One may choose HIGH first to get the same result.)
5. Select RETURN to get out of the menu.

MAIN MENU	ALARM LIMIT	ALARM SOUND : ON	
PREV MENU			

TEMPERATURE ALARM LIMIT			
RETURN	UNITS	LOW	HIGH
TEMP	°C	30.0	42.0

ALARM SOUND

The menu which decides activation of warning signs and message displays when the respiration alarm is on.

MAIN MENU	ALARM LIMIT	ALARM SOUND : ON	
PREV MENU			

MAIN MENU	ALARM LIMIT	ALARM SOUND : OFF	
PREV MENU			

UNIT SELECT

Able to select unit with °C, °F.

MAIN MENU	ALARM		UNIT SELECT: °C

MAIN MENU	ALARM		UNIT SELECT: °F

10. PRINT

10.1 Print

Printer and Heat Sensitivity Paper
Function and Setup Menu

10.2 Paper Change

10.1 Print

Printer and Heat Sensitivity Paper

A printer used to print data onto thermal paper.

Size of the thermal paper roll: 580mm wide x 380mm in diameter any thermal paper of same size can be used for the printer.

Side View of Printer



Function and Setup Menu

MAIN MENU	PRINTER SPEED: 25mm/s		WAVE FORM1: ECG
PREV MENU	WAVE FORM3: SPO2	WAVE FORM2: RESP	

1. Press the PRINT Key for continuous printing.

2. Select Printing Speed 25, 50 mm/s.

MAIN MENU	PRINTER SPEED: 25mm/s		WAVE FORM1: ECG
PREV MENU	WAVE FORM3: SPO2	WAVE FORM2: RESP	


MAIN MENU	PRINTER SPEED: 50mm/s		WAVE FORM1: ECG
PREV MENU	WAVE FORM3: SPO2	WAVE FORM2: RESP	

3. Set up ALARM PRINT in the MORE menu to activate ALARM during printing.



4. Data is printed in a selected wave form along with personal information of the animal.

3 channels select 3 parameters to print.

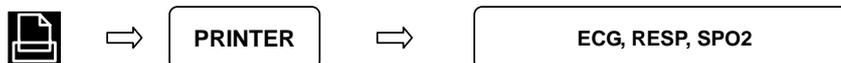

MAIN MENU	PRINTER SPEED: 50mm/s		WAVE FORM1: ECG
PREV MENU	WAVE FORM3: SPO2	WAVE FORM2: RESP	

MAIN MENU	PRINTER SPEED: 50mm/s	WAVE FORM1: ECG	> OFF ECG SPO2 RESP
PREV MENU	WAVE FORM3: SPO2		

MAIN MENU	PRINTER SPEED: 50mm/s		WAVE FORM1: ECG
PREV MENU	WAVE FORM3: SPO2	WAVE FORM2: RESP	

MAIN MENU	PRINTER SPEED: 50mm/s	WAVE FORM2: RESP	> OFF ECG SPO2 RESP
PREV MENU	WAVE FORM3: SPO2		

MAIN MENU	PRINTER SPEED: 50mm/s		WAVE FORM1: ECG
PREV MENU	WAVE FORM3: SPO2	WAVE FORM2: RESP	

MAIN MENU	WAVE FORM3: SPO2	OFF ECG SPO2 RESP	WAVE FORM1: ECG
PREV MENU		>	

10.2 Paper Change

1

Open the window of the printer.



2

Insert the paper roll offered with the product into the printing unit. Place the roll in a proper way so that the printed paper can roll out upwards.



3

Press the printer window until it is properly shut. Inaccurate shutting may cause failure in printing.



11. MESSAGE LIST

Function	Message	Details
ECG	LEAD FAULT	Cable is not properly connected.
SpO ₂	LEAD FAULT CHECK PROBE PULSE SEARCH POOR SIGNAL LOST PULSE ARTIFACT	Cable is not properly connected. Patient's finger is off the probe. Detection by the monitor of a pulse has ceased. The SpO ₂ signal is too low. The quality of the signal is questionable. The signal is patient's motion artifact
RESP	LEAD FAULT APNEA	Cable is not properly connected. APNEA gives an alarm.
NIBP	INFLATION FAILURE CHECK CUFF OVER PRESSURE DEFLATION FAILURE OVER TIME CUFF PRESSURE MEASUREMENT ERROR	Cuff hose is not properly connected. Cuff pressure is excessive. Cuff is bent, preventing deflation. Measure time exceeds the preset Level. Measure signal absent
TEMP	LEAD FAULT	Cable is not properly connected.
ALARM	ALARM VOL.OFF SILENCED ALARM PAUSE 5MIN	Alarm volume is off. Alarm key is pressed once Alarm key is pressed twice
TREND	NO PATIENT DATA	No patient's data input.
PRINT	NO PAPER	No paper in the printer
SETUP	BATTERY LOW	Low battery

12. DEFAULT SETTING VALUE

1. HORSE Mode

Alarm level

	High	Medium	Low	Message
Asystole	0			
Vfib/VTac	0			
V Tach	0			
HR		0		
NIBP		0		
SpO ₂			0	
SpO ₂ -Rate				0
RR				0
RR-Apnea				0
T(°C)				0

Parameter Limits

	Low	High
HR	60	150
NIBP-S	80	200
NIBP-M	50	170
NIBP-D	30	150
SpO ₂	90	100
SpO ₂ -Rate	60	150
RR(RESPI)	15	100
RR-Apnea	0	20
T(°C/ °F)	36.0/96.8	40.0/104.0
ST	-10.0	10.0
PVC	0	20

Display

Animal Age	0~2 years
Color format	Color
Primary ECG	II
Arrhythmia	Off
Detect Pace	Off
Print Waveform2	Off
Print Waveform3	Off
Alarm Print	Off
NIBP Auto	Off
NIBP Cuff Size	Large
RR(RESPI) Lead	II
Alarm Volume	Off
QRS Volume	Off
Pulse Volume	Off
ECG Lead Fault	Message
SpO ₂ Probe Off	Low Alarm
Units for Height	cm
Units for Weight	kg
Temperature Units	' F
NIBP Limit Type	Systolic
ECG Filter	Monitoring

2. DOG Mode

Alarm level

	High	Medium	Low	Message
Asystole	0			
Vfib/VTac	0			
V Tach	0			
HR		0		
NIBP		0		
SpO ₂			0	
SpO ₂ -Rate				0
RR				0
RR-Apnea				0
T(°C)				0

Parameter Limits

	Low	High
HR	60	160
NIBP-S	80	200
NIBP-M	50	170
NIBP-D	30	150
SpO ₂	90	100
SpO ₂ -Rate	60	160
RR(RESP)	15	100
RR-Apnea	0	20
T(°C/ F)	36.0/96.8	40.0/104.0
ST	-10.0	10.0
PVC	0	20

Display

Animal Age	0~2 years
Color format	Color
Primary ECG	II
Arrhythmia	Off
Detect Pace	Off
Print Waveform2	Off
Print Waveform3	Off
Alarm Print	Off
NIBP Auto	Off
NIBP Cuff Size	Medium
RR(RESP) Lead	II
Alarm Volume	Off
QRS Volume	Off
Pulse Volume	Off
ECG Lead Fault	Message
SpO₂ Probe Off	Low Alarm
Units for Height	cm
Units for Weight	kg
Temperature Units	° F
NIBP Limit Type	Systolic
ECG Filter	Monitoring

3. PUPPY Mode

Alarm level

	High	Medium	Low	Message
Asystole	0			
Vfib/VTac	0			
V Tach	0			
HR		0		
NIBP		0		
SpO ₂			0	
SpO ₂ -Rate				0
RR				0
RR-Apnea				0
T(° C)				0

Parameter Limits

	Low	High
HR	70	180
NIBP-S	80	200
NIBP-M	50	170
NIBP-D	30	150
SpO ₂	90	100
SpO ₂ -Rate	70	180
RR(RESP)	15	100
RR-Apnea	0	20
T(° C/ ° F)	36.5/97.7	39.5/103.1
ST	-10.0	10.0
PVC	0	20

Display

Animal Age	0~2 years
Color format	Color
Primary ECG	II
Arrhythmia	Off
Detect Pace	Off
Print Waveform2	Off
Print Waveform3	Off
Alarm Print	Off
NIBP Auto	Off
NIBP Cuff Size	Small
RR(RESP) Lead	II
Alarm Volume	Off
QRS Volume	Off
Pulse Volume	Off
ECG Lead Fault	Message
SpO₂ Probe Off	Low Alarm
Units for Height	cm
Units for Weight	kg
Temperature Units	° F
NIBP Limit Type	Systolic
ECG Filter	Monitoring

4. CAT Mode

Alarm level

	High	Medium	Low	Message
Asystole	0			
Vfib/VTac	0			
V Tach	0			
HR		0		
NIBP		0		
SpO ₂			0	
SpO ₂ -Rate				0
RR				0
RR-Apnea				0
T(° C)				0

Parameter Limits

	Low	High
HR	70	200
NIBP-S	80	200
NIBP-M	50	170
NIBP-D	30	150
SpO ₂	90	100
SpO ₂ -Rate	70	200
RR(RESP)	15	100
RR-Apnea	0	20
T(° C/ ° F)	36.5/97.7	39.5/103.1
ST	-10.0	10.0
PVC	0	20


Display

Animal Age	0~2 years
Color format	Color
Primary ECG	II
Arrhythmia	Off
Detect Pace	Off
Print Waveform2	Off
Print Waveform3	Off
Alarm Print	Off
NIBP Auto	Off
NIBP Cuff Size	Small
RR(RESPI) Lead	II
Alarm Volume	Off
QRS Volume	Off
Pulse Volume	Off
ECG Lead Fault	Message
SpO₂ Probe Off	Low Alarm
Units for Height	cm
Units for Weight	kg
Temperature Units	° F
NIBP Limit Type	Systolic
ECG Filter	Monitoring

13. TROUBLE SHOOTING


1. Noise in ECG

- Gel is dry
- Electrodes do not stick well to skin

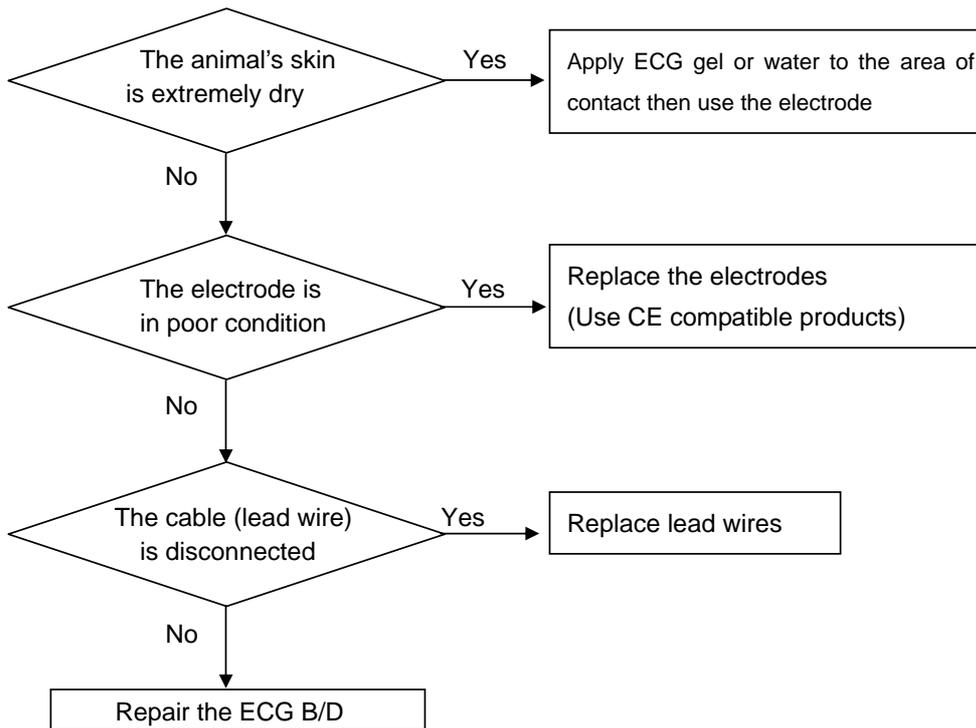


2. SpO₂ malfunction


Connectors are in bad condition?

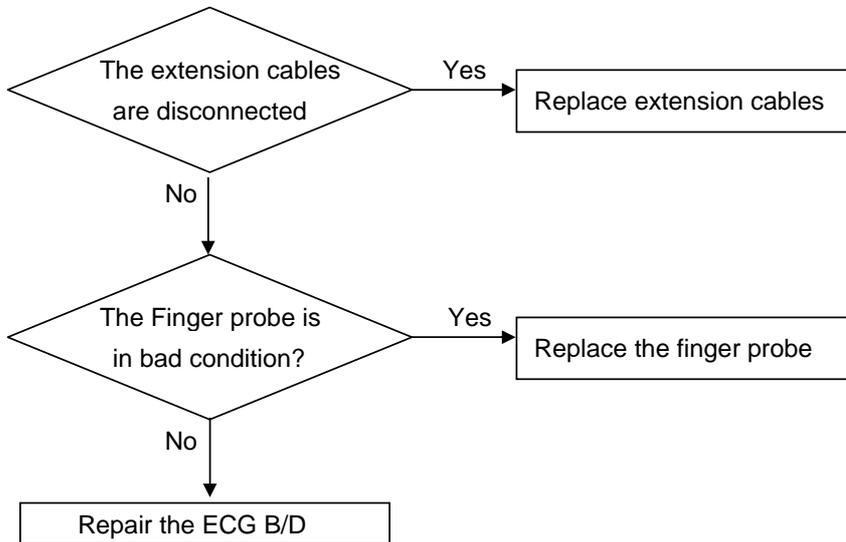


3. Temp malfunction

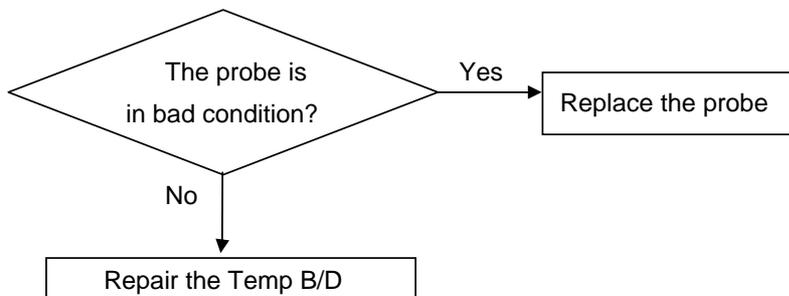

4. NIBP malfunction

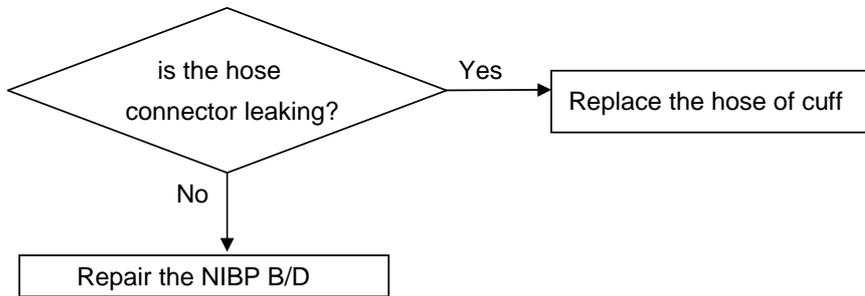


5. Abnormality in NIBP measurements

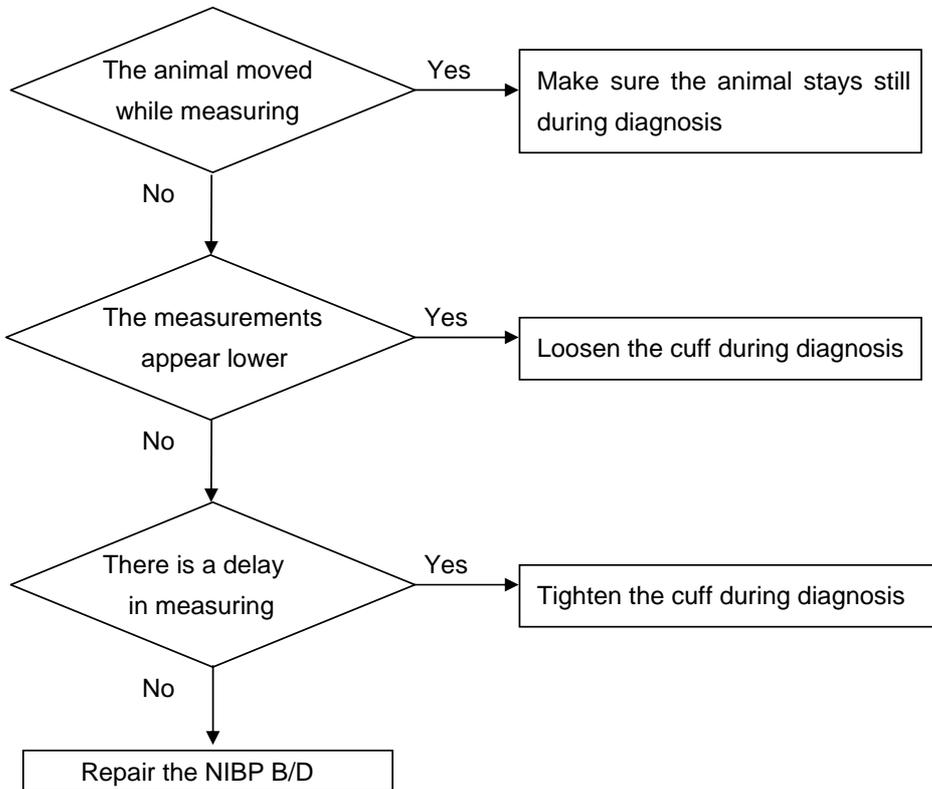


6. Failure in battery recharge


(the battery does not fully recharge in 6 hours or more)

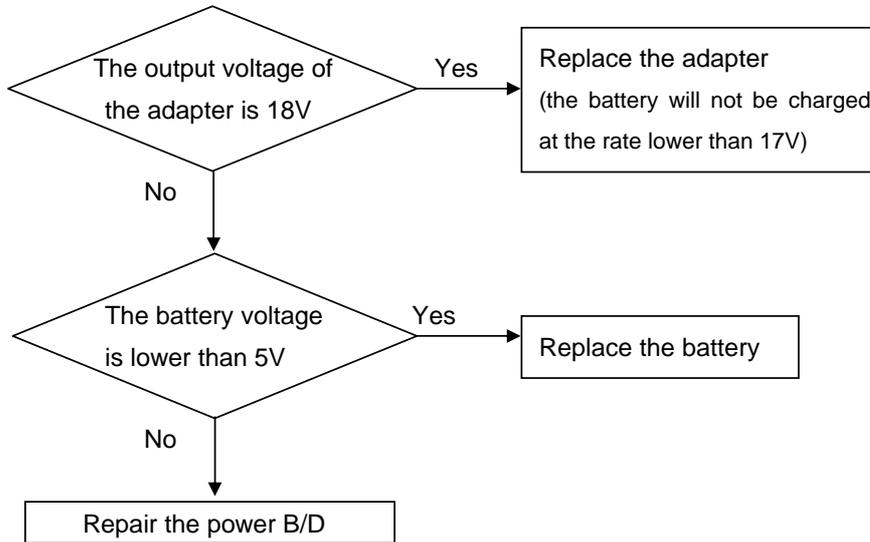


7. Power failure

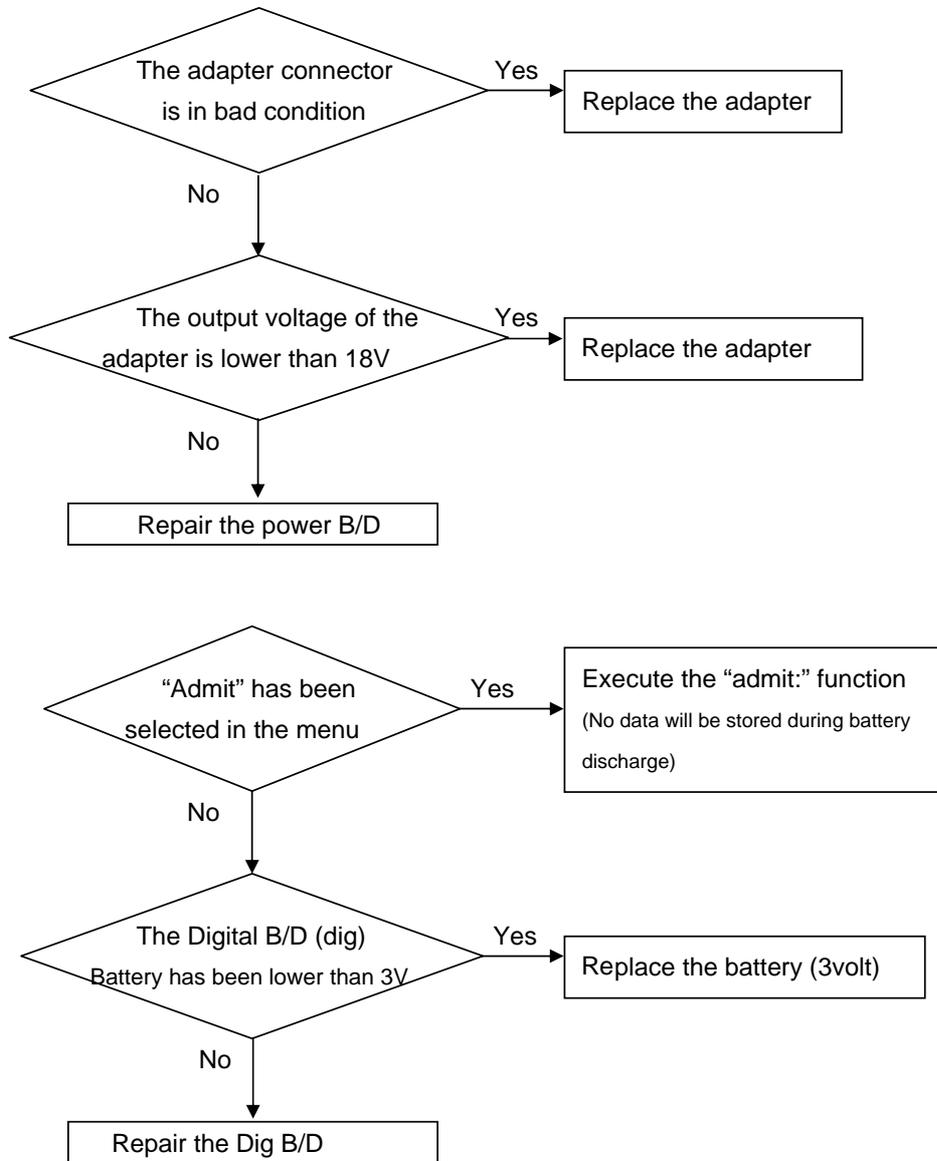

8. Periodic noises

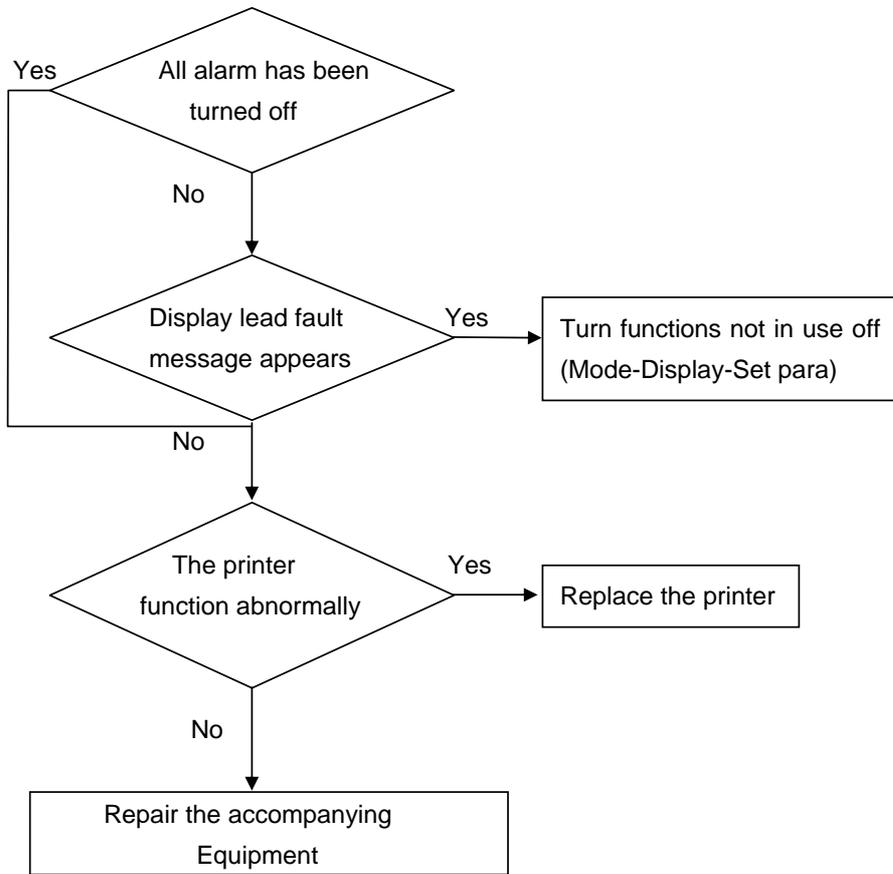

9. Print failure

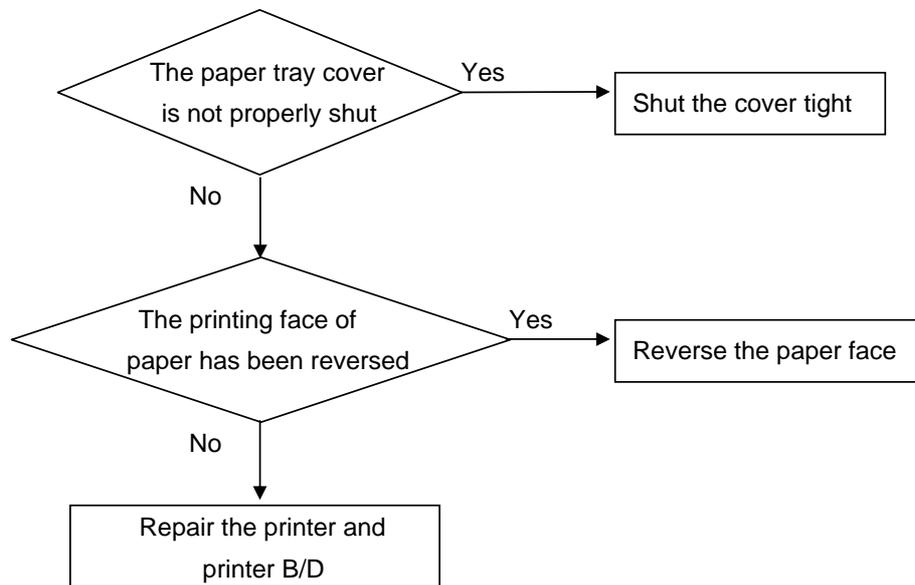

14. SPECIFICATION

Ease of use

Customization

Special Features

Monitor Environmental Specifications

Power adaptor

Monitor Performance Specifications

Graphical and Tabular Trends

SpO2 Performance Specifications

Respirations Performance Specifications

NIBP Performance Specifications

ECG Performance Specifications

Temperature Unit Performance Specifications

Accessories included

OPTION

Ease of use

- Battery operation
- Attached printer
- Table and graphic trend
- Nellcor SpO₂ sensor interchanges

Additional Function

- Able to use auto mobile power supply
- LAN Connection

Monitor Environmental Specifications

- Operating Temperature : 15°C to 30°C (59°F to 86°F)
- Storage Temperature : - 10°C to 60°C (14°F to 140°F)
- Humidity : 20% to 95% RH
- Operating Attitude : 70(700) to 106Kpa(1060mbar)

Power

- AC 100-240V (50/60Hz)
- Adapter 18 V, 2.5 A

Monitor Performance Specifications

- Screen : 7" TFT LCD (800×480)
- Indicators
 - Up to 3 wave patterns
 - 3 levels of alarm sound
 - Visual alarm
 - Pulse sound
 - handle flashing
 - Battery status
 - LED external power supply LED
- Interfaces
 - Vehicles power supply : 12 to 16 V DC, 3A max.
 - Generating power for LAN, Wireless LAN : 5.0V max 0.9A
- Battery

- Li-ion battery
- Battery status display
- Operating time : 2hours(with fully charged Battery)
- Thermal Printer(Optional) : internal printer
 - Speed : 25, 50 mm/sec
 - Paper width : 58 mm

Graphical and Tabular Trends

- Table Trend
 - Memory Storage : 128 hours
 - Data Interval : 1 minute
 - Display Interval : 1MIN, 5, 15, 30, 1HR
- Graphical Trend
 - Display Period : 30MINS, 60, 90, 3HRS, 6, 12

ECG capacity

- Lead : 3,5
- pulse rate range : 30 to 300 bpm
- pulse accuracy : ± 3 bpm
- Bandwidth : 0.5 Hz to 40 Hz
- Display Sweep Speed : 2.5mm / sec
- ECG size (Sensitivity) : 0.5, 1, 2, 4 mV/cm
- Lead-off Detection with display indicator
- Pace maker Detection Mode
- Differential Input Impedance : $> 5 \text{ M}\Omega$
- XCommon Mode Rejection Ratio : $> 90 \text{ dB}$ at 50 or 60 Hz
- DC Input Range : $\pm 5 \text{ mV}$
- Defibrillator Discharge : $< 5 \text{ s}$
- Defibrillation Artifact Recovery Time : $< 8 \text{ s}$

SpO₂ capacity

- Saturation Range : 0% to 100% oxygen proportion
- Pulse Rate Range : 30 to 254 bpm
- SpO₂ accuracy : 70% to 100% ± 2 digits, 0% to 69% unspecified

- pulse accuracy : ± 2 bpm
- Sensor Red 660nm, 2mW (typical)
Infrared 905nm, 2-2.4mW (typical)
- Minimum Signal: 0.05% modulation (Low perfusion level performance and
Amplitude limitation validation using FLUKE Index 2 Oximetry Simulator)

Respiration Performance Specifications

- Range : 5 to 120 breaths/min
- Accuracy : ± 3 breaths/min
- Display Sweep Speeds : 25mm/sec

NIBP capacity

- Technique : Oscillometric
- Measurement mode:
 - Manual : Single Measurement
 - Auto : automatic Intervals of 1MIN, 2, 3, 4, 5, 10, 15, 20, 30, 1HR, 2, 4, 8
- Pressure Display : 0 to 300 mmHg
- Blood Pressure Measurement Range: (LARGE, MEDIUM, and SMALL)
 - systolic : 60 to 250 mmHg
 - Mean Arterial Pressure : 40 to 235 mmHg
 - Diastolic : 30 to 220 mmHg

Temperature Unit Performance Specifications

- Range : 15°C to 45°C (59°F to 113°F)
- Accuracy : 25°C to 45°C $\pm 0.1^\circ\text{C}$, 15°C to 24°C $\pm 0.2^\circ\text{C}$
- Sensor : YSI 400 Series compatibility

Accessories Included:

- 3Lead animal cable (3WIRE-400) 1 EA
- Electrodes (ECG SENS - 400) 10 EA
- NIBP tubing, 3m long (NBPCBL-400) 1 EA
- Adult cuff, 25-35 Cm (ACUFF-430) 1 EA
- SpO2 extension cable 2m (SPCBL-400) 1 EA

- SpO₂ sensor (SPASENS – 430) 1 EA
- DC adapter, 18VDC, 2.5A (MW160 Made in AULT Co., Ltd.) 1 EA
- Printer Paper (PAPER-400) 2ROLL

Option

- Temperature sensor (skin) (TEMPSENS-430)
- 5 lead animal cable (5CBL-400, 5WIRE-430)

Abbreviations and Symbols

Abbreviations and symbols which you may encounter while reading this manual or using the monitor are listed below with their meanings.

Abbreviations

		A
A	amps	
AC	alternating current	
ARRHYTHM	arrhythmia	
ASYS	asystole	
Auto, AUTO	automatic	
AUX	Auxiliary	
aVF	left foot augmented lead	
aVL	left arm augmented lead	
aVR	right arm augmented lead	
		B
BPM	beats per minute	
		C
C	Celsius	
CAL	calibration	
cm, CM	centimeter	
		D
D	diastolic	
DC	direct current	
DEFIB, Defib	defibrillator	
DIA	diastolic	
		E
ECG	electrocardiograph	

EMC electromagnetic compatibility
 EMI electromagnetic interference
 ESU electrosurgical cautery unit

F

F Fahrenheit

G

g gram

H

HR heart rate, hour
 Hz hertz

I

ICU intensive care unit
 Inc incorporated

K

kg, KG kilogram
 kPa kilopascal

L

L liter, left
 LA left arm, left atrial
 LBS pounds
 LCD liquid crystal display
 LED light emitting diode
 LL left leg

M

M mean, minute
 m meter
 MEDI. Medium
 MIN, min minute
 MM, mm millimeters

MM/S millimeters per second
MMHG, mmHg millimeters of mercury
mV millivolt

N

NIBP noninvasive blood pressure

O

OR operating room

P

PVC premature ventricular complex

Q

QRS interval of ventricular depolarization

R

RA right arm, right atrial

RESP respiration

RL right leg

RR respiration rate

S

S systolic

sec second

SpO₂ arterial oxygen saturation from pulse oximetry

SYNC, Sync synchronization

SYS systolic

T

Temp, TEMP temperature

U**V**

V precordial lead

V	volt
V-Fib, VFIB	ventricular fibrillation
VTAC	ventricular tachycardia

W

X

X	multiplier when used with a number (2X)
---	---

Symbols

&	and
°	degree(s)
>	greater than
<	less than
–	minus
#	number
%	percent
±	plus or minus

Disposal

The product must not be disposed of along with other domestic waste. The users must dispose of this equipment by bringing it to a specific recycling point for electric and electronic equipment. For further information on recycling points contact the local authorities, the local recycling center or the shop where the product was purchased. If the equipment is not disposed of correctly, fines or penalties may be applied in accordance with the national legislation and regulations.

GIMA Warranty Conditions

Congratulations for purchasing a GIMA product. This product meets high qualitative standards both as regards the material and the production. The warranty is valid for 12 months from the date of supply of GIMA. During the period of validity of the warranty, GIMA will repair and/or replace free of charge all the defected parts due to production reasons. Labor costs and personnel traveling expenses and packaging not included. All components subject to wear are not included in the warranty. The repair or replacement performed during the warranty period shall not extend the warranty. The warranty is void in the following cases: repairs performed by unauthorized personnel or with non-original spare parts, defects caused by negligence or incorrect use. GIMA cannot be held responsible for malfunctioning on electronic devices or software due to outside agents such as: voltage changes, electro-magnetic fields, radio interferences, etc. The warranty is void if the above regulations are not observed and if the serial code (if available) has been removed, cancelled or changed. The defected products must be returned only to the dealer the product was purchased from. Products sent to GIMA will be rejected.



BIONET Co., Ltd. #1101, E&C Venture Dream Tower III, 38-21, Gi Guro-Gu, 152-719 Seoul - Korea



MGB Endoskopische Geräte GmbH Berlin, Schwarzschildstr. 6, 12489 Berlin - GERMANY